



Published in final edited form as:

*Nat Protoc.* 2006 ; 1(2): 729–748. doi:10.1038/nprot.2006.98.

## Chromatin immunoprecipitation and microarray-based analysis of protein location

Tong Ihn Lee<sup>1</sup>, Sarah E Johnstone<sup>1,2</sup>, and Richard A Young<sup>1,2</sup>

<sup>1</sup>Whitehead Institute for Biomedical Research, 9 Cambridge Center, Cambridge, Massachusetts 02142, USA.

<sup>2</sup>Department of Biology, Massachusetts Institute of Technology, Cambridge, Massachusetts 02139, USA.

### Abstract

Genome-wide location analysis, also known as ChIP-Chip, combines chromatin immunoprecipitation and DNA microarray analysis to identify protein-DNA interactions that occur in living cells. Protein-DNA interactions are captured *in vivo* by chemical crosslinking. Cell lysis, DNA fragmentation and immunoaffinity purification of the desired protein will co-purify DNA fragments that are associated with that protein. The enriched DNA population is then labeled, combined with a differentially labeled reference sample and applied to DNA microarrays to detect enriched signals. Various computational and bioinformatic approaches are then applied to normalize the enriched and reference channels, to connect signals to the portions of the genome that are represented on the DNA microarrays, to provide confidence metrics and to generate maps of protein-genome occupancy. Here, we describe the experimental protocols that we use from crosslinking of cells to hybridization of labeled material, together with insights into the aspects of these protocols that influence the results. These protocols require approximately 1 week to complete once sufficient numbers of cells have been obtained, and have been used to produce robust, high-quality ChIP-chip results in many different cell and tissue types.

### INTRODUCTION

Chromatin immunoprecipitation (ChIP) has become a widely used method for studying how proteins interact with the genome. The ChIP technique identifies physical interactions between proteins and DNA that occur within living cells<sup>1-3</sup> and ChIP can provide a whole-genome view of protein-DNA interactions when combined with DNA microarray analysis<sup>4-14</sup>. Our goal here is to share our ChIP-Chip protocol, suggest useful quality-control tests, identify the most likely sources of trouble and suggest solutions to problems with experiments. The protocol presented here has been in use for several years and has been used by many different individuals using a broad range of different mammalian cell types and tissues.

© 2006 Nature Publishing Group

Correspondence should be addressed to R.A.Y. (young@wi.mit.edu).

**COMPETING INTERESTS STATEMENTS** The authors declare competing financial interests (see the HTML version of this article for details).

Reprints and permissions information is available online at <http://npg.nature.com/reprintsandpermissions/>

## Overview of the ChIP-Chip technique

In a typical experiment (Figure 1), protein-DNA interactions are first captured *in vivo* by treating cells with a small, reversible crosslinker that is capable of rapidly diffusing into cells. After cell lysis, the DNA is fragmented and protein-DNA complexes can be enriched by immunoaffinity capture of the desired protein. As the DNA is fragmented prior to the capture, the DNA fragments that are most proximal to the binding event are enriched. Often, an aliquot of the input DNA is saved prior to immunoprecipitation for use as a reference sample. The crosslinks for both the enriched and input DNA are then reversed, and the DNA is purified away from RNA and protein components. The enrichment relative to a reference sample can then be measured by one of several techniques, such as PCR, comparative sequencing or DNA microarray analysis. For most DNA microarray experiments, the enriched and reference samples must be amplified first to provide sufficient material. In this protocol, the ends of the DNA are blunted and ligated to small DNA linkers that are subsequently used in priming PCR. This ligation-mediated PCR (LM-PCR) results in expanded pools of the enriched and reference DNA that are then differentially labeled and hybridized to the microarrays. The arrays are washed and then scanned. Various algorithms are then used to normalize the enriched and reference channels, to calculate the degree of enrichment and confidence metrics for each feature, to incorporate the information gained from having multiple probes representing contiguous regions of the genome and to generate maps of protein-genome occupancy.

## Applications of ChIP-Chip

The range of applications for ChIP-Chip is broad but still expanding. The method has recently been used to investigate DNA-binding transcription factors and the networks formed by these factors (reviewed in refs. 15–18), the distributions of chromatin-modifying machinery, histones, histone variants and histone modifications (reviewed in refs. 16–19), factors involved in DNA replication, DNA repair and DNA methylation (reviewed in refs. 20–22), and the linkage between transcription and nuclear-export machinery<sup>23–25</sup>. Immunoprecipitation has been coupled with microarrays to identify interactions between proteins and RNA<sup>26–30</sup> and ChIP has recently been used to study RNA-binding proteins that become associated with the genome through their interactions with RNA<sup>31</sup>. We anticipate that many new biological insights will emerge from protein-DNA and protein-RNA interactions that are identified using ChIP.

## Limitations of ChIP-Chip

There are several features of ChIP-Chip that currently limit its utility and make it a challenging technique<sup>32–34</sup>. The primary limitation is the quality of antibodies that are available for a protein of interest. The antibody, ideally, has high avidity and specificity for an epitope that is accessible in the crosslinked chromatin. Because this is not always the case, the investigator must be cautious about interpreting the lack of detection as reflecting the absence of an interaction. A second limitation is that a substantial number of cells ( $10^7$ – $10^8$ ) are generally needed to obtain a robust result. With current methods, using fewer cells can reduce signals, leading to an underestimation of protein-DNA interactions. High-quality, whole-genome DNA microarray platforms are expensive. As array manufacturers continue to make improvements in designs and manufacturing, more cost-effective experiments will be possible. Finally, the experimental process is long and can be challenging to troubleshoot.

## Experimental design

**Background**—The sections below describe several issues that we have found to be critical to consider in the early stages of planning and designing ChIP-Chip experiments.

**Antibodies**—ChIP-Chip is highly dependent on the antibody used for the immunoprecipitation. As individual antibodies can perform differently in the ChIP assay, it is helpful to screen a wide variety of available antibodies before launching a full-scale ChIP-Chip experiment. We generally start with antibodies, for which there is previous evidence for use in ChIP, or we use a variety of polyclonal and monoclonal antibodies from several sources.

The most reliable way to test the antibodies is to use them in small-scale ChIPs and to test for enrichment at a gene-specific level, although this approach depends entirely on the number of known binding sites, the strength of the evidence for the binding event and the comparability with your specific cell type, growth conditions and treatments. For small-scale ChIPs, we prepare material as described below, but instead of using all the material for a single reaction, we will split it into five to ten aliquots and use each aliquot in an immunoprecipitation. This lower amount of material seems adequate for gene-specific analysis, but does not usually work for microarray analysis. If no binding sites are known, it is much more difficult to screen for a useful antibody. Western blot analysis may be informative but is unreliable, as the antibody may recognize the denatured protein in a western blot but may still fail to recognize the native form *in vivo*.

If there are no antibodies for a specific protein, it is possible to use epitope-tagged versions of proteins, provided the tagged version does not radically affect protein function. We have used myc, TAP and FLAG epitope-tagged versions of various proteins with success.

**Reference samples**—Chromatin immunoprecipitation is an enrichment relative to a reference and not an absolute measurement, thus the choice of reference sample can be an important element of experimental design. Although it is most common to use unenriched, genomic DNA as the reference, it is sometimes more useful to hybridize an immunoprecipitation of one factor against an immunoprecipitation of a second factor to capture subtle effects. For instance, this approach can be useful when studying histone modifications. An immunoprecipitation against a core histone sub-unit is used as the reference to normalize the signal of the modification against the amount of histone. This can help to determine whether a perceived gain or loss of the modification is due to actual changes in modification or due to changes in histone density.

**Replicates**—ChIP-Chip, like many other forms of genome-wide analysis, requires the use of replicate samples to account for experimental noise that is inherent to high-throughput approaches and biological complexity due to the many factors influencing protein and nucleic-acid interactions. In particular, biological replicates (independently grown and treated pools of cells or tissue) or replicates with different antibodies against the same factor are valuable. The number of replicates needed depends on the noise that is inherent in the platform and the type of analysis to be performed, but triplicate experiments are preferred and duplicate experiments should be considered to be the minimum requirement whenever feasible.

**Experimental controls**—The use of a positive control and negative control for the experiment can be useful, both for trouble-shooting and for calibrating the analysis. The positive control should be a high-quality antibody against a protein that is likely to be present in any cell type. For example, in mammalian cells, we regularly use antibodies against the ubiquitous cell-cycle regulator E2F4. The most common negative control is an immunoprecipitation using a non-specific antibody that matches the isotype of the experimental antibody.

**Choice of arrays**—There are a number of options available for array hybridization, including self-printed oligonucleotide or PCR amplicon arrays, and commercial arrays and services. The choice of array platform ultimately depends on the need to balance performance, ease of use, resolution, cost of arrays and investment in equipment, and will vary from lab to lab. Due to the variety of options, a fully detailed protocol for our specific hybridization, washing, scanning and analysis steps is unlikely to be generally useful. However, we have included the relevant details for hybridization and washing and offer general guidelines that we try to follow for scanning and data analysis. This protocol was developed based on our work with both self-printed PCR amplicon arrays and commercially available oligonucleotide arrays, and reflects our most recent experience with oligonucleotide arrays<sup>13,14,35,36</sup>. A number of useful ChIP-Chip protocols have been published that describe protocols similar to this one, but which cover more specific aspects of dealing with different array platforms<sup>37–44</sup>.

**Normalization controls on the array**—As ChIP measures enrichment in one sample compared with another, it is essential to define an appropriate baseline for a particular experiment — one that reflects the signal from the background of all sequences from the whole genome and reflects lack of enrichment in the immunoprecipitation. This requirement can affect the sets of controls that are needed on the array. When there are a small number of features expected to show enrichment, the enrichment ratios of the majority of features can serve to represent the baseline level; the large size of this latter set strengthens the statistical significance of the enrichment. In other situations, for example, when examining RNA polymerase II binding at promoter regions using arrays focused only on promoter regions, we found that it was more difficult to determine the baseline enrichment that represented ‘no binding’. RNA polymerase II bound at a relatively high fraction of probes, making it more difficult to distinguish bound and unbound probes with as high a degree of statistical significance. In these cases, we found it useful to have a set of control probes designed against ‘gene deserts’, genomic regions that are devoid of known open reading frames and for which there is no evidence of noncoding transcripts. We assumed that these control probes would be generally unenriched for RNA polymerase II and could therefore serve as a normalization baseline. Similarly, we use control probes to subtract out a baseline level of non-specific hybridization signal and, for multislide arrays, control probes to normalize between slides.

**Data analysis and bioinformatics**—The ChIP-Chip analysis generates a large amount of data that requires some investment in statistical, computational and bioinformatics infrastructure.

Statistical analysis is required to account for experimental noise that is inherent to high-throughput approaches, but more importantly, to deal with the biological complexity of protein and nucleic-acid interactions. We expect that a protein of interest will be distributed across many different sites in a single cell in a manner that is dependent on protein concentration, its affinities for each site and other variables such as the presence of proteins that compete for similar sites. As a result, the output of ChIP-Chip is a distribution of enrichment ratios representing DNA or RNA occupancy, as opposed to two distinct sets of ‘bound’ and ‘not bound’. Statistical methods are needed to identify the statistically significant set of enrichment ratios, to provide a measure of confidence at the level of individual probe features and to incorporate data from multiple features to build models for binding that are both sensitive and specific. As there are a wide variety of options for statistical analysis<sup>7,12,45–51</sup> and the needs of individual labs will vary, a fully detailed description of our analysis methods is unlikely to be generally useful. However, we have included a basic outline of the approach that we have used in the past.

The volume and complexity of the data from ChIP-Chip analysis may require access to resources that are capable of creating or adapting programming code to manipulate and analyze the data. Even if most of the early steps are handled by a commercially available analysis package, there is almost always a need for additional analyses that need to be customized in some way. The capacity to handle and display the data using common programming languages, such as Visual Basic, MATLAB, Perl, R, Java or Python, is highly useful.

Along these lines, the need to be able to map ChIP-Chip data back to genomic location and genomic features often requires the development of bioinformatics resources, such as access to genome annotation and transcription-factor databases, and the ability to handle and cross-reference genome-wide expression data. Commonly used databases include the UCSC Genome Browser (<http://genome.ucsc.edu/>); the NCBI databases (<http://www.ncbi.nlm.nih.gov/>), particularly RefSeq, Gene and GEO datasets; the EMBL-EBI data-bases (<http://www.ebi.ac.uk/services/>); GO functional annotation (<http://www.geneontology.org/>); and TRANSFAC (<http://www.gene-regulation.com/pub/databases.html>).

### Protocol Tips

**Background**—The sections below describe a few key steps of the protocol for which we have some experience with variables that are particularly critical to the success of the experiment. The steps are presented in order, with those most likely to cause problems being presented first.

**Sonication**—Sonication is the most variable step in the process and will vary greatly depending on cell type, cell culture conditions used, quantity of cells, volume of sonication, degree of crosslinking and specifics of the sonicator being used. As a result, conditions must usually be optimized for each situation. Sonication fulfills two roles, solubilization of chromatin and shearing of chromatin, both of which are essential for successful ChIP analysis. Solubilization is achieved relatively easily. However, the degree of shearing requires more care. Undersonication results in a loss of resolution of binding events. Smaller DNA fragments allow for more precise localization of a specific binding event, as a smaller region of DNA will be pulled down in the immunoprecipitation. We have also found that oversonication must be avoided as it results in more noise in the microarray analysis and difficulty in identifying legitimate targets. As a guide, most protocols suggest sonication until reaching an average fragment size (different protocols suggest a value varying from 350 to 1,000 bp) or until the fragment sizes are within a certain range. Undersonication relative to these kinds of guidelines is easy to identify: the size of the reverse crosslinked, purified DNA fragments will be too large. Oversonication is more difficult to quantify. Most crosslinked material has a physical limit of sonication, and additional sonication past this point will not result in further visible shearing but will affect microarray results. Thus, it is possible that researchers trying to match previously published descriptions that were optimized for different conditions may inadvertently over-sonicate their DNA. As an alternative, we suggest sonication using the lowest power and time settings that result in sheared DNA, of which at least a quarter of the total DNA is sized from 200 to 600 bp in size, with less emphasis on the average size or overall spread of fragment sizes. One way to identify these conditions is to run a small time course of sonications. We will set up a sample of cells, as described below, with a default power setting but with an extended time of sonication. We then remove small aliquots (100  $\mu$ l out of 3.5 ml of material) at discrete time points. After crosslink reversal and purification, as described in the protocol below, the DNA can be run on a 2% agarose gel to estimate the degree of sonication. A sample sonication time course is shown in Figure 2.

**Immunoprecipitation conditions**—Antibody performance can be optimized by adjusting a number of variables, but we find that the simplest, most effective adjustment has been modifying the amount of salt in the immunoprecipitation. If gene-specific or microarray analysis indicates that antibody performance is poor (two-fold enrichment or less), but a positive control experiment indicates that the protocol is working well, it may be useful to perform a set of small-scale immunoprecipitations over a range of salt concentrations (typically starting with 0, 50, 100 and 250 mM sodium chloride). Prepare input DNA as described below but use only one-fifth to one-tenth of the material in an immunoprecipitation. This small-scale immunoprecipitation usually suffices for gene-specific analysis, but does not usually work for microarray analysis. If there is previous evidence indicating specific salt conditions that are optimal for a particular antibody, we will supplement the lysate with additional sodium chloride to approximate those conditions.

**Washing**—Following ChIP, it is possible to overwash the beads. This usually results in noise in the microarray analysis. In general, this occurs due to extensive wash times rather than too many washes. Five washes (as presented below) are adequate for most antibodies. If signal:noise ratios are low, increasing the number of washes (as high as 8) may help, although if you suspect low anti-body affinity, a lower number of washes (as few as 3) may be more appropriate.

**Crosslinking**—In practice, we have found that the crosslinking times and formaldehyde concentrations provided in this protocol are generally applicable and have been used successfully with a variety of proteins in human, mouse, yeast, fruit fly and zebra-fish model systems. But crosslinking times could theoretically be optimized for each protein. Insufficient crosslinking can result in the inability to capture protein-DNA contacts. Over-crosslinking can actually denature the protein of interest or cause over-aggregation, thus obscuring the epitope. Unfortunately, these problems are sometimes difficult to diagnose as the primary readout is the failure to detect enrichment, which could also result in problems at any number of other steps. One other sign would be unusual ease or difficulty in obtaining a range of sonicated material when testing a range of sonication conditions.

## MATERIALS

### REAGENTS

- 1 M HEPES-KOH, pKa 7.55 (Invitrogen/Gibco; Cat. no.: 15630-080)
- 5 M NaCl (Sigma S-5150)
- 500 mM EDTA (Invitrogen/Gibco; Cat. no: 15575-038)
- 500 mM EGTA (Sigma, Cat no: E3889)
- 37% formaldehyde (J.T. Baker; Cat. no: 2106-01) **! CAUTION** Formaldehyde is flammable; highly toxic by inhalation, contact with skin or if swallowed; causes burns; and is potentially carcinogenic. Formaldehyde should be used with appropriate safety measures, such as protective gloves, glasses and clothing, and adequate ventilation. Formaldehyde waste should be disposed of according to regulations for hazardous waste.
- 2.5 M glycine (Sigma; Cat. no: G8790)
- 10× Dulbecco's phosphate buffered saline (PBS) (Invitrogen/Gibco; Cat. no.: 14200-075)
- Dynabeads (Dyna) **▲ CRITICAL** Although it is possible to use other types of beads for the immunoprecipitation, we find that Dynabeads give us very low



background and are easy to use. We use Dynabeads that come pre-coated with secondary antibody (e.g., PanMouse IgG, which is a monoclonal human anti-mouse IgG) or pre-coated with protein G. The exact type of bead used will depend on the primary antibody being used)

- Block Solution: 1× PBS, 0.5% bovine serum albumin (BSA) (Sigma; Cat. no: A7906) ▲ **CRITICAL** Should be made fresh and kept cold.
- Primary antibody of choice
- 25× solution Complete Protease Inhibitor Cocktail: 1 tablet dissolved in 2 ml of double distilled water (can store up to 12 weeks at  $-20^{\circ}\text{C}$ ) (Roche; Cat. no.: 1 697 498)
- 50% glycerol (Sigma; Cat. no.: G5516)
- 10% Igepal CA-360 (Sigma; Cat. no.: I8896)
- 10% Triton X-100 (Sigma; Cat. no.: T8787)
- 1 M Tris-HCl, pH 8.0 (Invitrogen/Gibco; Cat. no.: 15568-025)
- 5% Na-Deoxycholate (Sigma; Cat. no.: D5760)
- 5% N-lauroylsarcosine (Fluka; Cat. no.: 61743) ▲ **CRITICAL** An ultrapure version of this reagent is not required for lysis steps, but is required for hybridizations. Lower grades can leave fluorescent residue on slides.
- 5 M LiCl (Sigma; Cat. no.: L7026)
- 10% SDS (Invitrogen/Gibco; Cat. no.: 15553-035)
- 10 mg ml<sup>-1</sup> RNaseA (Sigma; Cat. no.: R6513)
- Proteinase K solution (Invitrogen; Cat. no.: 25530-049)
- Phenol:chloroform:isoamyl alcohol (Fluka; Cat. no.: 77617) ▲ **CRITICAL** If this solution is old or is at low pH, there will be degradation of DNA. ! **CAUTION** Phenol is toxic when in contact with skin or if swallowed; causes burns; and is irritating to eyes, the respiratory system and skin. Chloroform is harmful by inhalation or if swallowed; is irritating to skin; and is potentially carcinogenic. Isoamyl alcohol is flammable; and is irritating to eyes, the respiratory system and skin. The phenol:chloroform:isoamyl alcohol should be used with appropriate safety measures, such as protective gloves, glasses and clothing, and adequate ventilation.  
The phenol:chloroform:isoamyl alcohol waste should be disposed of according to regulations for hazardous waste.
- 20 mg ml<sup>-1</sup> glycogen (Roche; Cat. no.: 901 393)
- 100% ethanol (Aaper)
- 80% ethanol (diluted from 100% ethanol)
- T4 DNA polymerase, 3 U μl<sup>-1</sup> (New England Biolabs; Cat. no.: M0203S)
- 10× NE Buffer 2 (supplied with T4 DNA Polymerase)
- 1 μg μl<sup>-1</sup> BSA (diluted from 10 mg ml<sup>-1</sup> stock supplied with T4 DNA Polymerase)
- 2.5 mM dNTP mix (2.5 mM each dNTP) (GE Healthcare Life Sciences; Cat. no.: 27-2035-03)

- 3 M sodium acetate (Sigma; Cat. no.: S-7899)
- T4 DNA ligase, 400 U  $\mu$ l (New England Biolabs; Cat. no.: M0202S)
- 5 $\times$  T4 DNA Ligase buffer (Invitrogen; Cat. no.: 46300-018)
- 40  $\mu$ M solution oligo oJW102 (MWG Biotech; 5'-GCGGTGACCCGGGAGATCTGAATTC-3')
- 40  $\mu$ M solution oligo oJW103 (MWG Biotech; 5'-GAATTCAGATC-3')
- 15  $\mu$ M Linker (see REAGENT SETUP)
- 10 $\times$  ThermoPol buffer (New England Biolabs; Cat. no.: B9004S)
- Double-distilled water
- AmpliTaq polymerase (Applied Biosystems)
- 7.5 M ammonium acetate (Sigma; Cat. no.: A2706)
- BioPrime Array CGH Genomic Labeling System (Invitrogen; Cat. no.: 18095-011; includes Labeling Module and Purification Module) ▲ **CRITICAL** We have experimented with other purification methods, but have found that these kits gave us the best yield while efficiently removing smaller DNA products. We have found that these small DNA products increased the background noise at lower signal intensity levels on our arrays.
- 1 mM Cy3-dUTP (PerkinElmer/NEN; Cat. no.: NEL578)
- 1 mM Cy5-dUTP (PerkinElmer/NEN; Cat. no.: NEL579)
- Formaldehyde Solution (see REAGENT SETUP)
- Lysis Buffer 1 (see REAGENT SETUP)
- Lysis Buffer 2 (see REAGENT SETUP)
- Lysis Buffer 3 (see REAGENT SETUP)
- Wash Buffer (RIPA) (see REAGENT SETUP)
- Elution Buffer: 50 mM Tris-HCl, pH 8.0, 10 mM EDTA, 1.0% SDS
- TE: 10 mM Tris-HCl, pH 8.0, 1 mM EDTA
- TE + NaCl: 10 mM Tris-HCl, pH 8.0, 1 mM EDTA, 50 mM NaCl (keep cold)
- 2 $\times$  Blunting Mix (See Step 44) ▲ **CRITICAL** Should be made fresh before use.
- 2 $\times$  Ligase Mix (See Step 52) ▲ **CRITICAL** Should be made fresh before use. Linkers used in this mix should be thawed on ice, used immediately and not re-used.
- LMPCR Mix A (See Step 58)
- LMPCR Mix B (See Step 58)
- Label Mix (Nucleotide Mix and Klenow are from Invitrogen labeling kit) (See Step 70)
- The following reagents are specific to our hybridization and washing protocol; as different laboratories will select different array platforms, these are intended as a reference: 200 ng  $\mu$ l<sup>-1</sup> sheared herring sperm DNA (Promega; Cat. no.: D1815); 8  $\mu$ g  $\mu$ l<sup>-1</sup> yeast tRNA (Invitrogen; Cat. no.: 15401-029); 1  $\mu$ g  $\mu$ l<sup>-1</sup> Cot-1 DNA (Invitrogen; Cat. no.: 15279-011; use species-specific Cot-1 DNA)



- 10× Control Targets (Agilent; Cat. no.: 5185-5976-P) (Reconstitute per manufacturer's directions. Used in protocol at 0.1×, not 1× as recommended by manufacturer)
- 500 mM Na-MES, pH 6.9 (Sigma; Cat. no.: M5287) (Adjust pH with 10 N sodium hydroxide) ▲ **CRITICAL** MES can go off. As a precaution, make small batches of the stock solution and use within 3 months. Store in the dark at 4 °C.
- 100% formamide (Sigma; Cat. no.: F9037) ! **CAUTION** Formamide is harmful by inhalation, in contact with skin or if swallowed; and causes burns. Formamide should be used with appropriate safety measures, such as protective gloves, glasses and clothing, and adequate ventilation. The formamide waste should be disposed of according to regulations.
- 100% acetonitrile (J.T. Baker; Cat. no.: 9017-03) ! **CAUTION** Acetonitrile is highly flammable; and is toxic by inhalation, in contact with skin or if swallowed. Acetonitrile should be used with appropriate safety measures, such as protective gloves, glasses and clothing, and adequate ventilation. The acetonitrile waste should be disposed of according to regulations.
- Hybridization Control Mix (See Step 80)
- Hybridization Buffer Master Mix (See Step 82) ▲ **CRITICAL** Should be made fresh before use.
- 20× SSPE (Invitrogen; Cat. no.: 15591-027)
- Wash I: 6× SSPE, 0.005% *N*-lauroylsarcosine
- Wash II: 6× SSPE
- Wash III (Stabilization and Drying Buffer; Agilent; Cat. no.: 5185-5979) ! **CAUTION** This solution contains acetonitrile, which is flammable and hazardous and should be used in a fume hood.

## EQUIPMENT

- Cell strainer, 100 µm nylon (if using tissues) (e.g., BD Falcon)
- Rotator (e.g., Fisher Hematology/Chemistry Mixer)
- Magnetic particle collector (MPC; Dynal)
- Sonicator (e.g., Misonix 3000 sonicator equipped with microtip)
- Phase lock, heavy, pre-dispensed into 2.0 ml microfuge tubes (Eppendorf) (Although these could be replaced with standard liquid organic solvent extractions, the ease of use factor makes these ideal)
- NanoDrop ND-1000 spectrophotometer (This is highly useful as it allows you to assay DNA concentrations and dye incorporations by using low volumes (1.5 µl) of sample)
- Thermal cycler (e.g., MJ Research PTC-225) (For large-scale set-ups, it is useful to have a cycler that can handle 96-well plates)
- Swinging bucket centrifuge, variable temperature (e.g., Sorvall Legend RT)
- Rocking platform (e.g., Bellco Rocker Platform)
- Liquid nitrogen and appropriate container

- 65 degree oven (e.g., Techne Hybridiser HD-1B) (Ovens or warm rooms are preferred for longer incubations as the evenly applied heat limits condensation forming on tube lids)

## REAGENT SETUP

**15  $\mu$ M Linker**—Mix 250  $\mu$ l of 1 M Tris-HCl (pH 8.0), 375  $\mu$ l of 40  $\mu$ M oJW102 and 375  $\mu$ l of 40  $\mu$ M oJW103. Make 50 or 100  $\mu$ l aliquots in PCR tubes. Place tubes in thermal cycler. Set up and run the following program: 95 °C for 5 min; 70 °C for 1 min; ramp down to 4 °C (0.4 °C/min); 4 °C Hold. Remove linkers and store at -20 °C. **▲ CRITICAL** Linkers must be thawed on ice and used immediately. The linkers are labile and will dissociate, causing low-efficiency ligation. Aliquots of partially used linkers should not be re-frozen and re-used.

**Formaldehyde Solution**—Should be made fresh. Consists of 50 mM HEPES-KOH, pH 7.5, 100 mM NaCl, 1 mM EDTA, 0.5 mM EGTA and 11% formaldehyde.

**Lysis Buffer 1**—Add protease inhibitors just before use, filter and keep cold. Consists of 50 mM HEPES-KOH, pH 7.5, 140 mM NaCl, 1 mM EDTA, 10% glycerol, 0.5% NP-40, 0.25% Triton X-100, 1 $\times$  protease inhibitors

**Lysis Buffer 2**—Add protease inhibitors just before use, filter and keep cold. Consists of 10 mM Tris-HCl, pH 8.0, 200 mM NaCl, 1 mM EDTA, 0.5 mM EGTA, 1 $\times$  protease inhibitors

**Lysis Buffer 3**—Add protease inhibitors just before use, filter and keep cold. Consists of 10 mM Tris-HCl, pH 8.0, 100 mM NaCl, 1 mM EDTA, 0.5 mM EGTA, 0.1% Na-Deoxycholate, 0.5% *N*-lauroylsarcosine, 1 $\times$  protease inhibitors

**Wash Buffer (RIPA)**—Keep cold. Consists of 50 mM HEPES-KOH, pKa 7.55, 500 mM LiCl, 1 mM EDTA, 1.0% NP-40, 0.7% Na-Deoxycholate

## EQUIPMENT SETUP

**Sonicator**—We most commonly use a Misonix 3000 sonicator with microtip. The settings we typically begin with are power setting 7 (~35 W) and 12 cycles (where each cycle is a 30 s burst of sonication, followed by a 90 s pause). The range of settings will vary based on cell type, cell number, growth conditions and crosslinking. For comparison, the range of settings we have used for specific experiments is: power settings from 6 to 9; cycle number from 8 to 15; burst times of 20 to 30 s and pause times of 60 to 90 s.

**Array scanning**—As the specifics of scanning depend on the platform chosen and the scanner available, it is less useful to provide a detailed protocol on how to scan arrays using any one specific set-up. However, there are a few guidelines that we have found to be useful for two-color hybridizations. In selecting scanner settings, the primary variable we change is photomultiplier tube (PMT) gain, which adjusts the sensitivity of the PMT that converts photons to a digitized electrical signal for the Cy3 and Cy5 channels.

A higher PMT gain creates a brighter image. We try to select PMT settings for each of the lasers in a two-color scanner that maximize signal intensity from features, without creating excessive background signal from the glass surface. Our current array platform has extremely low background signal, so our main concern is maximizing signal intensity. PMT gains are also set such that the distributions of signal intensities for each channel are similar. In addition, the settings should not result in a large number of signals that exceed, and are

therefore thresholded at, the maximum detection range of the scanner. This can result in loss of detection of enrichment for features with high signal intensities.

## PROCEDURE

### Formaldehyde crosslinking cells ● TIMING Day 1, 1–4 h

1| For formaldehyde-crosslinking of cells, three procedures are described that can be used, depending on whether your cells are: adherent to culture apparatus (A); growing in suspension (B); or obtained fresh from tissues (C). The procedure given in (C) has been used with mouse and human tissues. Specifics vary from tissue to tissue, and this extraction procedure will therefore require some optimization — these steps are, therefore, general guidelines for extracting and crosslinking cells from tissues.

! **CAUTION** See REAGENTS for precautions when using formaldehyde.

#### ? TROUBLESHOOTING

##### A. For adherent cells:

- i. Use  $5 \times 10^7 - 1 \times 10^8$  cells for each immunoprecipitation.
- ii. Add 1/10 volume of fresh 11% Formaldehyde Solution to plates. Swirl briefly.
- iii. Incubate cells with Formaldehyde Solution for 10 min at room temperature.
- iv. Add 1/20 volume of 2.5 M glycine to quench formaldehyde.
- v. Rinse cells twice with 10 ml of 1× PBS. Harvest cells using silicon scraper.
- vi. Pool cells in 50 ml conical tubes and spin at 700g for 5 min at 4 °C. Discard supernatant and resuspend pellet in 10 ml of 1× PBS per  $10^8$  cells with gentle inversion (cells may stick to pipettes at this stage).
- vii. Aliquot  $5 \times 10^7 - 1 \times 10^8$  cells to individual 15 ml conical tubes and spin at 700g for 5 min at 4 °C. Discard supernatants.

##### B. For suspension cells:

- i. Use  $5 \times 10^7 - 1 \times 10^8$  cells for each immunoprecipitation.
- ii. Add 1/10 volume of fresh 11% Formaldehyde Solution to flasks. Swirl briefly.
- iii. Incubate cells with Formaldehyde Solution for 20 min at room temperature.
- iv. Add 1/20 volume of 2.5 M glycine to quench formaldehyde.
- v. Pool cells in required number of 50 ml conical tubes and spin at 700g for 5 min at 4 °C. Discard supernatant.
- vi. Resuspend cells in 50 ml of 1× PBS with gentle inversion (cells may stick to pipettes at this stage). Spin at 700g for 5 min at 4 °C to pellet cells. Discard supernatant. Repeat once. Resuspend final cell pellet in 10 ml of 1× PBS per  $10^8$  cells.

- vii. Aliquot  $5 \times 10^7 - 1 \times 10^8$  cells to individual 15 ml conical tubes and spin at 700g for 5 min at 4 °C. Discard supernatants.

**C. For tissues:**

- i. Harvest tissues. For human tissues, this involves an arrangement to obtain transplant-grade tissue that has been released for research and will depend on the tissue. Whenever possible, tissues are crosslinked as close to collection as possible and then shipped. If not possible, tissues are shipped on ice, allowed to recover briefly and then crosslinked and processed as described below. For mouse tissues, this entails euthanizing the mice and then immediately harvesting the desired tissues.
- ii. Mince tissues very finely. Transfer to 50 ml tube and add 2× volume of 1× PBS. Add 1/10 volume of fresh 11% Formaldehyde Solution.
- iii. Incubate cells with Formaldehyde Solution for 10 min at room temperature. Swirl tubes occasionally. Add 1/20 volume of 2.5 M glycine to quench formaldehyde. Immediately place cells on ice.
- iv. Homogenize cells. We typically use douncing or hand-held, mechanical homogenizers. Pass material through a 100-µm nylon cell strainer.
- v. Pool cells in required number of 50 ml conical tubes and spin at 1,100g for 5 min at 4 °C. Discard supernatant.
- vi. Resuspend cells in 50 ml 1× PBS with gentle inversion (cells may stick to pipettes at this stage). Spin at 1,100g for 5 min at 4 °C to pellet cells. Discard supernatant. Resuspend final cell pellet in 10 ml of 1× PBS per  $10^8$  cells.
- vii. Aliquot  $5 \times 10^7 - 1 \times 10^8$  cells to individual 15 ml conical tubes and spin at 1,350g for 5 min at 4 °C. Discard supernatants.

2| Flash-freeze cells in liquid nitrogen and store pellets at  $-80$  °C.

■ **PAUSE POINT** Once cells are crosslinked, they may be stored frozen at  $-80$  °C indefinitely.

**Preparing magnetic beads ● TIMING Day 2, 7–8 h**

3| Add 100 µl of Dynal beads to a 1.5 ml microfuge tube. Set up 1 tube of beads per immunoprecipitation. Add 1 ml Block Solution.

▲ **CRITICAL STEP** Steps 3–9 are all performed at 4 °C.

4| Collect beads using Dynal MPC. Place tubes in rack. Allow beads to collect on side of tube. This should take approximately 15 s. Invert rack once or twice to help collect all beads. Remove supernatant with aspirator or pipettor.

5| Add 1.5 ml Block Solution and gently resuspend beads. This can be done by removing the tubes from the rack and inverting the tubes 10–20 times until the beads are evenly resuspended or by removing the magnetic strip from the rack and inverting the rack, with the tubes still in place — either 10–20 times or until the beads are evenly resuspended. Collect beads as above (Step 4). Remove supernatant with aspirator or pipettor.

- 6| Wash beads in 1.5 ml Block Solution, as in Step 5, one more time.
- 7| Resuspend beads in 250  $\mu$ l Block Solution and add 10  $\mu$ g of antibody. Incubate at 4  $^{\circ}$ C for a minimum of 6 h, or overnight, on a rotator.

■ **PAUSE POINT** This step can be extended to overnight at 4  $^{\circ}$ C. If so, beads should be prepared starting on Day 1.

? **TROUBLESHOOTING**

- 8| Wash beads three times in 1 ml Block Solution, as described in Step 5.
- 9| Resuspend each aliquot of beads in 100  $\mu$ l Block Solution.

**Cell sonication** ● **TIMING Day 2, 2 h**

- 10| Remove frozen cell pellets from  $-80^{\circ}$ C and resuspend each pellet of  $\sim 10^8$  cells in 10 ml of Lysis Buffer 1. Rock at 4  $^{\circ}$ C on platform rocker for 10 min.
- 11| Spin at 1,350g for 5 min at 4  $^{\circ}$ C in a table-top centrifuge. Discard supernatant.
- 12| Resuspend each pellet in 10 ml of Lysis Buffer 2. Rock gently at room temperature for 10 min.
- 13| Pellet nuclei in table-top centrifuge by spinning at 1,350g for 5 min at 4  $^{\circ}$ C. Discard supernatant.
- 14| Resuspend each pellet in each tube in 3.5 ml Lysis Buffer 3. Transfer cells to tubes for sonication. We currently prefer to use the bottom half of a standard polypropylene, 15 ml conical tube for sonication. We cut the tube into two pieces at the 7 ml mark and discard the upper half. Tubes can be covered with parafilm while setting up. You may also wish to save a 50  $\mu$ l aliquot for use as a pre-sonication control.
- 15| Using a clamp, position tube so the sonicator probe sits approximately 0.5–1.0 cm above the bottom of the tube. Take care that the probe is centered and does not contact the sides of the tube.

▲ **CRITICAL STEP** Probe positioning can affect whether the solution foams or not during sonication. Typically, foaming indicates that the sonicated DNA will be poorly sheared.

? **TROUBLESHOOTING**

- 16| Immerse tube in an ice-water bath. This is most easily done by keeping the sonicator probe and tube fixed, placing the bath on a height-adjustable platform and raising it into position.
- 17| Sonicate suspension. Samples should be kept in an ice-water bath during sonication. To decrease foaming, initially set output power to 4 and increase manually to final power during first burst.

▲ **CRITICAL STEP** If there is significant foaming, we recommend removing all bubbles by centrifugation at 20,000g followed by gentle resuspension of all material, leaving no foam bubbles.

? **TROUBLESHOOTING**

- 18| Add 1/10 volume of 10% Triton X-100 to sonicated lysate. Split into two 1.5 ml centrifuge tubes. Spin at 20,000g for 10 min at 4  $^{\circ}$ C to pellet debris.
- 19| Combine supernatants from the two 1.5 ml centrifuge tubes in a new 15 ml conical tube for immunoprecipitation.

▲ **CRITICAL STEP** From here on, keep lysates on ice.

20| Save 50  $\mu$ l of cell lysate from each sample as input DNA. Store at  $-20^{\circ}\text{C}$ . At least one input DNA aliquot should be kept per batch of sonicated lysate. Note that the effects of the sonication and the resulting distribution of fragment sizes can only be checked at Step 42 after crosslink reversal and purification of DNA.

■ **PAUSE POINT** Lysates can be frozen at  $-80^{\circ}\text{C}$  and used the next day.

### Chromatin immunoprecipitation ● TIMING Day 2, 15 min and overnight incubation

21| Add 100  $\mu$ l antibody/magnetic bead mix from Step 9 to cell lysates from Step 19. Incubate overnight on rotator at  $4^{\circ}\text{C}$ .

▲ **CRITICAL STEP** Antibody performance can be affected by the amount of salt in the reaction. If there is previous evidence indicating specific salt conditions that are optimal for a particular antibody, we will supplement the lysate with additional sodium chloride to approximate those conditions. In this case, the salt concentrations of wash buffers used in Steps 25–26 should also be adjusted.

#### ? TROUBLESHOOTING

### Wash ● TIMING Day 3, 2 h

22| Set up 1.5 ml microfuge tubes on ice, one tube for each immunoprecipitation. Transfer half of each immunoprecipitation from the 15 ml conical tube to its assigned microfuge tube.

▲ **CRITICAL STEP** All wash steps are performed at  $4^{\circ}\text{C}$ . Buffers should be kept cold.

23| Collect beads using Dynal MPC. Place tubes in rack. Allow beads to collect on side of tube. This should take approximately 15 s. Invert rack once or twice to help collect all beads. Remove supernatant with aspirator or pipettor, changing tips between samples.

24| Add second half of immunoprecipitation to tube. Let tubes sit again in magnetic stand to collect the beads. Remove supernatant with aspirator or pipettor, changing tips between samples.

25| Add 1 ml Wash Buffer (RIPA) to each tube and gently resuspend beads. This can be done by removing the tubes from the rack and inverting the tubes 10–20 times or by removing the magnetic strip from the rack and inverting the rack, with tubes still in place — 10–20 times or until the beads are evenly resuspended. Collect beads. Remove supernatant by aspirator or pipettor. Repeat this wash four more times, changing tips between washes.

#### ? TROUBLESHOOTING

26| Wash once with 1 ml TE + 50 mM NaCl.

27| Spin at 960g for 3 min at  $4^{\circ}\text{C}$  and remove any residual TE buffer.

### Elution ● TIMING Day 3, 30–45 min

28| Add 210  $\mu$ l of Elution Buffer and elute material from beads by incubating tubes in a  $65^{\circ}\text{C}$  water bath for 15 min. Vortex briefly every 2 min. This incubation can be extended as long as 30 min, which can help improve recovery of the eluate.

29| Spin down beads at 16,000g for 1 min at room temperature.



**30|** Remove 200  $\mu\text{l}$  of supernatant and transfer to new tube.

■ **PAUSE POINT** Material can be frozen at  $-20\text{ }^{\circ}\text{C}$  and stored overnight.

### Crosslink reversal ● **TIMING Day 3, 6 h or overnight**

**31|** Reverse crosslink the immunoprecipitation DNA from Step 30 by incubating at  $65\text{ }^{\circ}\text{C}$  for a minimum of 6 h. This incubation is usually done in an oven so that the tube is heated evenly and there is less condensation formed.

#### ? **TROUBLESHOOTING**

**32|** Thaw 50  $\mu\text{l}$  of input DNA reserved after sonication (Step 20), add 150  $\mu\text{l}$  (3 volumes) of elution buffer and mix. Reverse crosslink this input DNA by incubating at  $65\text{ }^{\circ}\text{C}$  as in Step 31. From this point, every tube of immunoprecipitation or input DNA is considered to be a separate tube or sample for later processing steps.

■ **PAUSE POINT** If the wash was performed early in the day, the crosslink reversal can be stopped after a minimum of 6 h and the material can be frozen at  $-20\text{ }^{\circ}\text{C}$  and stored overnight. If the wash was performed later in the day, the crosslink reversal can be extended up to 15 h and performed overnight with little effect on the reaction.

▲ **CRITICAL STEP** Longer times of crosslink reversal (18 h or more) usually result in increased noise in the microarray analysis.

#### ? **TROUBLESHOOTING**

### Purification of DNA ● **TIMING Day 4, 6 h**

**33|** Add 200  $\mu\text{l}$  of TE to each tube to dilute SDS in Elution Buffer. We have found that high levels of SDS can inhibit RNase activity.

**34|** Add 8  $\mu\text{l}$  of 10  $\text{mg ml}^{-1}$  RNaseA (0.2  $\text{mg ml}^{-1}$  final concentration), mix by inverting the tube several times and incubate at  $37\text{ }^{\circ}\text{C}$  for 2 h.

■ **PAUSE POINT** Material can be frozen at  $-20\text{ }^{\circ}\text{C}$  and stored overnight.

**35|** Add 4  $\mu\text{l}$  of 20  $\text{mg ml}^{-1}$  Proteinase K (0.2  $\mu\text{g ml}^{-1}$  final concentration) and mix by inverting the tube several times and incubate at  $55\text{ }^{\circ}\text{C}$  for 2 h.

**36|** Add 400  $\mu\text{l}$  phenol:chloroform:isoamyl alcohol (P:C:IA), vortex and separate phases with 2 ml Heavy Phaselock tube (follow instructions provided by Eppendorf).

▲ **CRITICAL STEP** If the P:C:IA solution is old or is at low pH, there will be degradation of DNA, causing noise in the microarray analysis and loss of detection of valid targets.

! **CAUTION** see REAGENTS for precautions when using phenol and chloroform.

#### ? **TROUBLESHOOTING**

**37|** Transfer aqueous layer to new centrifuge tube containing 16  $\mu\text{l}$  of 5M NaCl (200 mM final concentration) and 1.5  $\mu\text{l}$  of 20  $\mu\text{g } \mu\text{l}^{-1}$  glycogen (30  $\mu\text{g}$  total).

**38|** Add 800  $\mu\text{l}$  EtOH. Incubate for 30 min at  $-20\text{ }^{\circ}\text{C}$ .

**39|** Spin at 20,000g for 10 min at  $4\text{ }^{\circ}\text{C}$  to pellet DNA. Wash pellets by adding 500  $\mu\text{l}$  of 80% EtOH, vortexing to resuspend pellet and spinning again at 20,000g for 5 min at  $4\text{ }^{\circ}\text{C}$ .

**40|** Remove any remaining 80% EtOH. Spin the tubes briefly to collect any remaining liquid and remove liquid with a pipettman, avoiding the pellet. Let tubes air dry until

pellets are just dry: pellets should still retain a moist appearance. Resuspend each pellet in 70  $\mu\text{l}$  of 10 mM Tris-HCl, pH 8.0.

▲ **CRITICAL STEP** Overdrying of these pellets can make them difficult to resuspend, or liable to flake and peel away from the side of the tube.

#### ? TROUBLESHOOTING

41| Save 15  $\mu\text{l}$  of immunoprecipitation sample for future use. This material can be used to perform gene-specific PCR confirmation of microarray results.

42| Measure DNA concentration of input DNA with NanoDrop and dilute input DNA to 100 ng  $\mu\text{l}^{-1}$ . Run purified input DNA on a 2% agarose gel to check distribution of fragment sizes and digestion of RNA. The immunoprecipitation DNA will usually be too dilute to measure or visualize at this point. If possible, this is a good opportunity to check for enrichment in the immunoprecipitation sample using gene-specific PCR.

■ **PAUSE POINT** Material can be frozen at  $-80\text{ }^{\circ}\text{C}$  and stored indefinitely.

#### ? TROUBLESHOOTING

### T4 polymerase blunting of ends ● TIMING Day 4, 2 h

43| Set up 1.5 ml microfuge tubes for each immunoprecipitation and input DNA sample. For each input DNA sample, mix 2  $\mu\text{l}$  (200 ng) input DNA and 53  $\mu\text{l}$  ddH<sub>2</sub>O in an individual tube. For each immunoprecipitation, aliquot 55  $\mu\text{l}$  into an individual tube.

▲ **CRITICAL STEP** Steps 43–45 should be performed on ice, unless otherwise stated.

44| Make 2× Blunting Master Mix according to the table below:

Component	Amount (per reaction)	Final amount/concentration
10× NEB buffer 2	11.0 $\mu\text{l}$	1×
1 $\mu\text{g } \mu\text{l}^{-1}$ BSA	5.5 $\mu\text{l}$	50 $\mu\text{g ml}^{-1}$
2.5 mM dNTP mix	4.4 $\mu\text{l}$	100 $\mu\text{M}$ each dNTP
T4 DNA polymerase (3 U $\mu\text{l}^{-1}$ )	0.5 $\mu\text{l}$	1.5 U
ddH <sub>2</sub> O	33.6 $\mu\text{l}$	
TOTAL volume	55.0 $\mu\text{l}$	

45| Add 55  $\mu\text{l}$  of 2× Blunting Master Mix to all tubes. Incubate for 20 min at 12  $^{\circ}\text{C}$  in a water bath that is set up in a 4  $^{\circ}\text{C}$  cold room.

▲ **CRITICAL STEP** The reaction should be kept at 12  $^{\circ}\text{C}$  and promptly removed after incubation. Theoretically, the low temperature, short reaction time and high concentrations of nucleotides favor the polymerase activity of the T4 DNA polymerase and not the exonuclease activity, resulting in blunt ends and not 5' overhangs. If no reliable water bath/cold room set-up exists, set up PCR tubes at Step 43 and use a thermal cycler to incubate reactions at 12  $^{\circ}\text{C}$ . Use of PCR tubes will require transferring the samples to microfuge tubes for Steps 46–47.

#### ? TROUBLESHOOTING

46| Add 11.5  $\mu\text{l}$  of 3 M sodium acetate and 0.5  $\mu\text{l}$  of 20  $\mu\text{g } \mu\text{l}$  glycogen (10  $\mu\text{g}$  total).

47| Add 120  $\mu\text{l}$  P:C:IA, vortex samples and transfer samples to room temperature Phaselock tubes. Extract once (follow instructions provided by Eppendorf).

▲ **CRITICAL STEP** If the P:C:IA solution is old or is at low pH, there will be degradation of DNA, causing noise in the microarray analysis and loss of detection of valid targets.

! **CAUTION** See REAGENTS for precautions when using phenol and chloroform.

#### ? TROUBLESHOOTING

48| Transfer aqueous layer to new centrifuge tube containing 250  $\mu$ l EtOH. Incubate for 30 min at  $-20^{\circ}\text{C}$ .

49| Spin at 20,000g for 10 min at  $4^{\circ}\text{C}$  to pellet DNA. Wash pellets by adding 500  $\mu$ l of 80% EtOH, vortexing to resuspend pellet and spinning again at 20,000g for 5 min at  $4^{\circ}\text{C}$ .

50| Remove any remaining 80% EtOH. Spin the tubes briefly to collect any remaining liquid and remove liquid with a pipetteman, avoiding the pellet. Let tubes air dry until pellets are just dry: pellets should still retain a moist appearance.

▲ **CRITICAL STEP** Overdrying of these pellets can make them difficult to resuspend, or liable to flake and peel away from the side of the tube.

#### ? TROUBLESHOOTING

51| Carefully resuspend each pellet in 25  $\mu$ l  $\text{H}_2\text{O}$  with pipetting. Adding the water and then allowing the pellets to sit on ice for some time can help. Chill on ice after resuspending.

■ **PAUSE POINT** Material can be frozen at  $-80^{\circ}\text{C}$  and stored indefinitely.

### Blunt-end ligation ● **TIMING Day 4, 1 h and overnight incubation; Day 5, 2 h**

52| Make 2 $\times$  Ligase Master Mix according to the table below:

Component	Amount (per reaction)	Final amount/concentration
5 $\times$ T4 DNA ligase buffer	10.0 $\mu$ l	1 $\times$
15 $\mu$ M linker	6.7 $\mu$ l	2 $\mu$ M
T4 DNA ligase (400 U $\mu\text{l}^{-1}$ )	0.5 $\mu$ l	200 U
dd $\text{H}_2\text{O}$	7.8 $\mu$ l	
TOTAL volume	25.0 $\mu$ l	

▲ **CRITICAL STEP** Steps 52–57 should be done on ice, unless otherwise stated.

53| Add 25  $\mu$ l Ligase Master Mix to 25  $\mu$ l of sample and incubate for 16 h in a  $16^{\circ}\text{C}$  water bath.

54| Add 6  $\mu$ l of 3 M sodium acetate and 130  $\mu$ l EtOH. Incubate for 30 min at  $-20^{\circ}\text{C}$ .

55| Spin at 20,000g for 10 min at  $4^{\circ}\text{C}$  to pellet DNA. Wash pellets by adding 500  $\mu$ l of 80% EtOH, vortexing to resuspend pellet and spinning again at 20,000g for 5 min at  $4^{\circ}\text{C}$ .

56| Remove any remaining 80% EtOH. Spin the tubes briefly to collect any remaining liquid and remove liquid with a pipetteman, avoiding the pellet. Pellet may appear as a translucent smear on the side of the tube. Let tubes air dry until pellets are just dry: pellets should still retain a moist appearance.

▲ **CRITICAL STEP** Overdrying of these pellets can make them difficult to resuspend, or liable to flake and peel away from the side of the tube.

? **TROUBLESHOOTING**

57| Carefully resuspend each pellet in 25  $\mu\text{l}$   $\text{H}_2\text{O}$  with pipetting. Adding the water and then allowing the pellets to sit on ice for some time can help. Store on ice after resuspending.

■ **PAUSE POINT** Material can be frozen at  $-80^\circ\text{C}$  and stored indefinitely.

**Ligation-mediated PCR** ● **TIMING Day 5, 3–6 h**

58| Make LMPCR Mix A and LMPCR Mix B according to the tables below. Mixes should be prepared on ice.

LMPCR Mix A		
Component	Amount (per reaction)	Final amount/concentration
10 $\times$ ThermoPol buffer	4.0 $\mu\text{l}$	1 $\times$
2.5 mM dNTP mix	5.0 $\mu\text{l}$	250 $\mu\text{M}$ each dNTP
40 $\mu\text{M}$ oJW102 oligonucleotide	1.25 $\mu\text{l}$	1 $\mu\text{M}$
ddH <sub>2</sub> O	4.75 $\mu\text{l}$	
TOTAL volume	15.0 $\mu\text{l}$	
LMPCR Mix B		
Component	Amount (per reaction)	Final amount/concentration
10 $\times$ ThermoPol buffer	1.0 $\mu\text{l}$	1 $\times$
Taq polymerase (5 U $\mu\text{l}^{-1}$ )	0.5 $\mu\text{l}$	2.5 U
ddH <sub>2</sub> O	8.5 $\mu\text{l}$	
TOTAL volume	10.0 $\mu\text{l}$	

59| Add 15  $\mu\text{l}$  of LMPCR Mix A to each 25  $\mu\text{l}$  sample, on ice. Transfer tubes to the thermal cycler.

60| If an additional expansion step is required to increase the amount of material for labeling, follow (A) below. If no expansion is required, run the program described in (B).

**A. Expansion required**

- i. Program the thermocycler as follows:

Cycle	Denature	Anneal	Extend
1	55 $^\circ\text{C}$ , 4 min	72 $^\circ\text{C}$ , 3 min	
2	95 $^\circ\text{C}$ , 2 min		
3–15	95 $^\circ\text{C}$ , 30 s	60 $^\circ\text{C}$ , 30 s	72 $^\circ\text{C}$ , 1 min

- ii. Start program. Midway through the 55 °C step of cycle 1, add 10  $\mu$ l Mix B to each tube to hot-start reactions. If necessary, pause the program so the tubes remain at 55 °C while adding Mix B.
- iii. After the 15-cycle program has finished, remove the samples and add 475  $\mu$ l of double-distilled water. This diluted material can then be used as template for subsequent LMPCR reactions (Steps 58–60 (B)). A total of 5  $\mu$ l of diluted material should be adequate for a new 50- $\mu$ l reaction. Adjust the amount of water in the recipes listed above accordingly.

▲ **CRITICAL STEP** It is strongly recommended that you do not use this expansion until comfortable with the procedure. It is important to first calibrate the number of cycles in the first-round PCR to ensure that 15 cycles of amplification are still in the linear range. We used both real-time PCR and array hybridization of material expanded once with 25 cycles versus material expanded with the 15/25 protocol to ensure that our results were not affected by the first-round expansion.

#### B. No expansion required

- i. If no expansion step is required, program the thermocycler as follows:

Cycle	Denature	Anneal	Extend
1	55 °C, 4 min	72 °C, 3 min	
2	95 °C, 2 min		
3–27	95 °C, 30 s	60 °C, 30 s	72 °C, 1 min
28	72 °C, 5 min		

The sample can be kept at 4 °C by adding a HOLD step at the end of the program.

- ii. Start program. Midway through the 55 °C step of cycle 1, add 10  $\mu$ l of Mix B to each tube to hot-start reactions. If necessary, pause the program so the tubes remain at 55 °C while adding Mix B.

**61** After PCR is completed, transfer reactions to individual 1.5 ml Eppendorf tubes. Pool samples if appropriate. For each 50  $\mu$ l of PCR reaction, add 25  $\mu$ l of 7.5 M ammonium acetate and 225  $\mu$ l of cold EtOH. Incubate for 30 min at –20 °C.

**62** Spin at 20,000g for 10 min at 4 °C to pellet DNA. Wash pellets with 500  $\mu$ l of 80% EtOH. Wash pellets by adding 500  $\mu$ l of 80% EtOH, vortexing to resuspend pellet and spinning again at 20,000g for 5 min at 4 °C.

**63** Dry pellets as described above (Step 56) and resuspend each pellet in 50  $\mu$ l H<sub>2</sub>O.

#### ? TROUBLESHOOTING

**64** Measure DNA concentration with NanoDrop (use 10-fold dilutions, if necessary) and normalize all samples to 500 ng  $\mu$ l<sup>-1</sup>. Run a 2% agarose gel to check size of LMPCR products.

■ **PAUSE POINT** Material can be frozen at –80 °C and stored indefinitely.

**Cy3/Cy5 labeling of IP/WCE material ● TIMING Day 6, 6 h**

**65|** Preheat thermal cycler to 95 °C.

**66|** For each labeling reaction, aliquot 2 µl (1 µg) of LMPCR-amplified immunoprecipitated or input DNA into a PCR tube.

**67|** Add 35 µl of 2.5× random primer solution and 38 µl of double-distilled water. This is a random-primed, Klenow-based extension protocol derived from Invitrogen's CGH kit. Our protocol varies from the instructions provided by Invitrogen in both reaction volume and reagent concentrations. We perform 20 reactions per '30 reaction' Invitrogen kit. A pair of reactions (one for immunoprecipitation and one for its reference input DNA) yields enough material for 1–2 hybridizations. To scale-up for more arrays, we increase number — not volume — of individual reactions.

**68|** Vortex samples for 30 s.

**69|** Place tubes in pre-heated thermal cycler. After 5 min at 95 °C, immediately transfer samples to an ice-water slurry bath for 10 min to flash-cool reagents.

**70|** While reagents are cooling, make up Label Mix. One mix is required for immunoprecipitation samples (usually Cy5). A second mix is required for input DNA samples (usually Cy3). If reagents permit, a dye swap can provide a useful comparison (a second set of replicates where the immunoprecipitate is labeled with Cy3 and the input DNA is labeled with Cy5).

**Label Mix:**

Component	Amount (per reaction)	Final amount/concentration
10× dUTP Nucleotide mix	8.2 µl	1×
Cy5 (or Cy3)	1.5 µl	17 µM
Klenow (40 U µl <sup>-1</sup> )	1.5 µl	60 U
ddH <sub>2</sub> O	1.8 µl	
TOTAL volume	13.0 µl	

**71|** Vortex samples for 30 s.

**72|** Add 13 µl of Label Mix to each sample. Pipette up and down multiple times to mix reagents.

**73|** Incubate samples for 3 h at 37 °C.

**▲ CRITICAL STEP** Keep samples in dark as Cy dyes are sensitive to exposure to light. Other incubation times and conditions have been tried, but we find that these conditions offer the best yield, while limiting the accumulation of smaller DNA fragments (<100 bp) that can contribute to noise in the microarray analysis.

**74|** Add 9 µl stop buffer to each well and mix by pipetting gently.

**75|** Clean up DNA using the DNA purification columns provided with the CGH kit. Follow manufacturer's instructions.

**? TROUBLESHOOTING**

**76|** After eluting labeled DNA from the column into a microfuge tube, spin tubes at 15,000g for 2 min at room temperature to pellet any fine particulates. Carefully transfer supernatants to fresh tubes. Repeat if necessary to ensure that the solution is clear.



▲ **CRITICAL STEP** It is particularly important to thoroughly remove the particulates at this step. Particulates from the purification column can bind to some slide surfaces, making it difficult, if not impossible, to successfully scan the slide.

#### ? TROUBLESHOOTING

77| Measure DNA concentration and dye incorporation with NanoDrop. Optional: Run a 2% agarose gel to check size of labeling products.

■ **PAUSE POINT** If this DNA is going to be used within the next day, it can be stored at  $-20^{\circ}\text{C}$ . If the DNA will be stored for longer, we recommend storing it as a pellet. Precipitate the DNA by adding 25  $\mu\text{l}$  of ammonium acetate and 300  $\mu\text{l}$  of ice-cold EtOH and placing the tubes at  $-80^{\circ}\text{C}$  for 10 min. Spin at 20,000g for 10 min at  $4^{\circ}\text{C}$  to pellet DNA. Wash pellets with 500  $\mu\text{l}$  of 80% EtOH. Remove excess ethanol. Allow pellets to air-dry briefly and store tubes at  $-80^{\circ}\text{C}$  in a light-proof box or wrapped in aluminum foil.

#### ? TROUBLESHOOTING

### Array hybridization ● TIMING Day 6, 1–2 h, plus 40-h incubation

78| For each hybridization, combine 5  $\mu\text{g}$  of the Cy5-labeled immunoprecipitated DNA and 5  $\mu\text{g}$  of the appropriate Cy3-labeled reference DNA. Keep all tubes on ice. The hybridization can be performed with less DNA (we have used as little as 3  $\mu\text{g}$  of each Cy5-labeled and Cy3-labeled), but better signal:noise ratios are seen with larger amounts of DNA. If working with frozen pellets, resuspend each pellet first in 50  $\mu\text{l}$  of double-distilled water and then combine however much you need for the hybridization.

79| For each hybridization, bring the combined DNA up to a total volume of 120  $\mu\text{l}$ . Keep all tubes on ice. If the DNA solutions are dilute, this volume can be extended as high as 170  $\mu\text{l}$  and water can be removed from the hybridization buffer mix as necessary.

80| Make Hybridization Control Master Mix.

#### Hybridization Control Master Mix

Component	Amount (per reaction)	Final amount/concentration
Sheared herring sperm DNA (200 ng $\mu\text{l}^{-1}$ )	4.0 $\mu\text{l}$	1.6 $\mu\text{g ml}^{-1}$
Yeast tRNA (8 $\mu\text{g ml}^{-1}$ )	5.0 $\mu\text{l}$	80 $\mu\text{g ml}^{-1}$
Cot-1 DNA (1.6 $\mu\text{g ml}^{-1}$ )	10.0 $\mu\text{l}$	20 $\mu\text{g ml}^{-1}$
Control oligonucleotides (1 $\times$ )	5.0 $\mu\text{l}$	0.1 $\times$
TOTAL volume	24.0 $\mu\text{l}$	

81| Add 24  $\mu\text{l}$  of Hybridization Control Master Mix to each hybridization. Mix by pipetting gently. Keep all tubes on ice.

#### Hybridization Buffer Master Mix

Component	Amount (per reaction)	Final amount/concentration
ddH <sub>2</sub> O	50.0 $\mu\text{l}$	

**Hybridization Buffer Master Mix**

Component	Amount (per reaction)	Final amount/concentration
500 mM Na-MES, pH 6.9	50.0 $\mu$ l	50 mM
5 M NaCl	50.0 $\mu$ l	500 mM
500 mM EDTA	6.0 $\mu$ l	6 mM
5 % <i>N</i> -lauroylsarcosine	50.0 $\mu$ l	0.5%
Formamide	150.0 $\mu$ l	30%
TOTAL volume	356.0 $\mu$ l	

**82** Make Hybridization Buffer Master Mix. This solution should be made fresh before use and kept at room temperature.

**!** CAUTION See REAGENTS for precautions when working with formamide.

**83** Add 356  $\mu$ l of Hybridization Buffer Master Mix to each hybridization. Mix by pipetting gently. After mixing, do not return samples to ice. Instead, leave tubes at room temperature until all samples have been processed. Keep all samples protected from light.

**84** Heat samples for 3 min in a heat block that has been set to 95 °C. Cover with foil to prevent exposure to light.

**85** Transfer tubes to a 40 °C heat block for 15 min. Cover with foil to prevent exposure to light. At this point, start preparing hybridization chambers for assembly, aiming to handle 4–6 slides at a time. If there are more slides than this, samples can sit at this 40 °C step for at least 1 h.

**86** After 15 min, spin tubes at 14,000g for 45 s to pellet any particulates.

**? TROUBLESHOOTING**

**87** Assemble hybridizations according to manufacturer's instructions. Use 490  $\mu$ l of hybridization mix per hybridization.

**88** Incubate at 40 °C in rotating oven for 40 h. Ensure free rotation of liquid throughout the chamber. We find the best results are obtained with an incubation time of 38–42 h.

**Array washing ● TIMING Day 8, 1–2 h**

**89** Before disassembling hybridization chambers, prepare Wash I (6 $\times$  SSPE, 0.005% *H*-lauroylsarcosine) and Wash II (6 $\times$  SSPE). Filter Wash I and Wash II before use. Place Wash III at 37 °C to dissolve solute. Swirling the bottle of Wash III occasionally will help to dissolve the solute. Allow the solution to cool to room temperature before use.

**!** CAUTION As Wash III is volatile, avoid temperatures that are higher than 45 °C.

**? TROUBLESHOOTING**

**90** Prepare and wash all glassware and equipment that is required for slide disassembly (follow manufacturer's instructions) and for washes described in Steps 91–95. We usually use standard glass 10-slide microscope slide racks and dishes.

**? TROUBLESHOOTING**

**91** After hybridization, disassemble hybridization chambers individually following the manufacturer's instructions. Batches of 4–6 slides are convenient. Slides are placed in

slide racks in dishes containing a magnetic stir bar (thin enough to spin freely when slides and slide rack are added) and sufficient Wash I to cover the slides.

#### ? TROUBLESHOOTING

92| Once all slides for this batch are in the dish, place dish on stir plate and gently stir for 5 min at room temperature. On a typical stir plate with settings from 1 to 10, use setting 3 or 4. The surface of the liquid should barely dimple. During this wash, prepare two dishes with 270 ml of Wash II, one with a stir bar.

#### ? TROUBLESHOOTING

93| After the 5-min wash in Wash I, remove rack from Wash I, quickly dip it into a dish containing Wash II and immediately transfer to the second dish containing 270 ml of Wash II and stir bar.

#### ? TROUBLESHOOTING

94| Place this dish on a stir plate and gently stir for 5 min at room temperature. Use stir plate settings as described in Step 92. During this wash, prepare the washes for the next step in the fume hood. Pour 270 ml of 100% acetonitrile (FPLC grade) into a small rectangular dish. Pour 250 ml (one half bottle) of room temperature Wash III into a second, small rectangular dish, add a stir bar and place on a stir plate. Cover both dishes when not in use to prevent evaporation.

! CAUTION See REAGENTS for precautions when using acetonitrile.

#### ? TROUBLESHOOTING

95| Turn off stir plate and carry dish containing the slides to the hood. Turn on stir plate in the hood to begin agitating Wash III. Remove rack, quickly dip it into dish containing 100% acetonitrile, and immediately transfer to a final dish containing Wash III and a rotating stir bar.

#### ? TROUBLESHOOTING

96| After 15 s, turn off the stir plate. This allows particulates to settle, and helps prevent them from sticking to your slide as you pull the slides out. Wait another 15 s.

97| Slowly and smoothly remove rack and slides from Wash III — it should take approximately 10 s to lift the rack from the solution (the slides should be dry at this point). If you move the rack abruptly as you are removing it, you may create visible lines of Wash III precipitate. If this happens, re-dip the slide momentarily and try to pull it out again.

#### ? TROUBLESHOOTING

98| Scan immediately. As fluorescent dyes are sensitive to exposure to atmosphere, it is best to scan slides as soon as possible after washing. Various ozone-scavenging treatments, such as Wash III, extend the life of hybridized slides, but still eventually fail. Slides can be stored in vacuum-sealed bags or a high nitrogen environment for short times (2–5 h) until needed.

#### ? TROUBLESHOOTING

### Array analysis

99| As described above, each lab is likely to have specific preferences for analysis of the array. A very brief, general outline is provided here as an example and a more detailed explanation is available in Lee *et al.*<sup>14</sup>. Scan slides (we use an Agilent DNA microarray scanner) at wavelengths of 532 nm for input (Cy3) and 635 nm for enriched (Cy5).

**100|** Use GenePix feature-identification software to automatically identify features and to extract data for fluorescence at each wavelength by calculating the average intensities for all pixels within the feature boundaries.

**101|** Subtract local background signal for each feature. Normalize the signal intensities of the slide across the entire set of slides using signals from a set of common features that are printed on every slide in the set. Then, subtract the signal intensities from a set of negative controls for hybridization.

**102|** Median normalize the data from the two channels and generate a log ratio of signal intensities for each feature.

**103|** Process the log ratios using a whole chip error model<sup>52</sup>, which incorporates signal intensity and background noise to calculate a *P*-value for the statistical significance of the enrichment ratio that is seen at each feature. We also calculate the *P*-value for enrichment for sets of three adjacent features and incorporate information about the enrichment at both individual features and sets of features to identify genomic regions that qualify as being bound by the protein in question.

### ? TROUBLESHOOTING

**104|** Identify distinct binding events by collapsing adjacent bound regions.

### ● TIMING

Day 1: Steps 1–2 (approximately 1–4 h, depending on numbers of cells or tissues being collected)

Day 2: Steps 3–9 (approximately 8 h, should be started first)

Day 2: Steps 10–21 (approximately 2 h and overnight incubation)

Day 3: Steps 22–32 (approximately 8–9 h, or approximately 2 h and overnight incubation)

Day 4: Steps 33–53 (approximately 9 h and overnight incubation)

Day 5: Steps 54–64 (approximately 6–7 h)

Day 6: Steps 65–88 (approximately 7–8 h and 40-h hybridization)

Day 8: Steps 89–98 (approximately 1–2 h and scanning time)

## TROUBLESHOOTING

Troubleshooting advice can be found in Table 1.

## ANTICIPATED RESULTS

### Step 42 (i)

Measure DNA concentration of input DNA (from Step 20) with NanoDrop. Input DNA samples should be approximately 70  $\mu\text{l}$  total volume with a DNA concentration of 300–500  $\text{ng } \mu\text{l}^{-1}$ . Immunoprecipitation samples are usually too dilute to detect reliably and measurements are sometimes skewed due to the glycogen that is added to the sample, but we have pooled and precipitated several immunoprecipitates and, after adjusting for changes in volume, expect a DNA concentration of 2–6  $\text{ng } \mu\text{l}^{-1}$ . The immunoprecipitated DNA concentration may be higher if the immunoprecipitation is directed against an abundant protein.

**Step 42 (ii)**

Run a 2% agarose gel of the input DNA that is saved prior to immunoprecipitation (Step 20). This allows you to visually check the distribution of fragment sizes that result from sonication, and to confirm digestion of RNA. Examples of different distributions of fragment sizes resulting from different degrees of sonication are shown in Figure 2. We would first try ChIP with the fragment distribution that is seen with 12 cycles of sonication, as this is the lowest amount of sonication that results in a majority of fragments, ranging from 200–600 bp. The fragment distribution that is seen with eight cycles of sonication is also likely to yield reasonable ChIP results, but this material has generally larger fragments that could reduce the ability to resolve binding events. The fragment distribution that is seen with 16 cycles of sonication is similar to the 12-cycle material, but we would default to the lower amount of sonication. Note that different cell types, growth conditions and crosslinking conditions will change the appearance of similar sonication time courses. However, the basic principle of finding a setting that minimizes sonication while maximizing fragments in the range of 200–600 bp should still be a useful starting point.

It can also be useful to calculate the fragment-size distribution more precisely with laser densitometry or a device such as the Agilent 2100 Bioanalyzer. These values are useful for computational algorithms that combine knowledge of the size distribution with genomic probe locations, probe characteristics and enrichment ratios to develop more sophisticated models of protein-DNA binding.

**Step 42 (iii)**

Whenever possible, you should check that the immunoprecipitation was successful by looking for enrichment of known positive controls in the immunoprecipitated DNA sample. This is most easily done by performing PCR using primers for known positive control and negative control regions and running the products on an agarose gel. Compare the ratio of DNA for the positive control versus the negative control for both the immunoprecipitated DNA (use 2  $\mu$ l) and the input DNA (use dilution series). The immunoprecipitation positive:negative ratio should be at least twice the input DNA positive:negative ratio. It is important to perform a serial dilution of the input DNA PCR reactions to ensure that all of the reactions are in the linear range of detection. Ideally, we run the reactions as a multiplex PCR so the expected positive and the negative control are being tested in the same PCR reaction (adjust primer concentrations so the input DNA positive:negative ratio is close to 1). This is not always possible with every combination of primers.

**Step 64 (i)**

Measure DNA concentration with NanoDrop. DNA samples should be approximately 50  $\mu$ l total volume with a DNA concentration of 100–350  $\text{ng } \mu\text{l}^{-1}$ . All samples, immunoprecipitated and input DNA should be comparable in DNA amounts.

**Step 64 (ii)**

Run a 2% agarose gel of the LMPCR DNA. DNA products should be centered around 200–400 bp. All samples, immunoprecipitated and input DNA should be comparable in distribution of fragment sizes and DNA amounts. An example of LMPCR DNA run on an agarose gel is shown in Figure 3.

**Step 77 (i)**

Measure DNA concentration and dye incorporation with NanoDrop. DNA samples should be approximately 50  $\mu$ l total volume with a DNA concentration of 100–150  $\text{ng } \mu\text{l}^{-1}$  and 4–10  $\text{pmol } \mu\text{l}^{-1}$  of incorporated dye.

**Step 77 (ii)**

Run a 2% agarose gel of the labeled, purified DNA. DNA products should be centered around 200–300 bp. There should not be significant amounts of DNA fragments that are smaller than 100 bp in the purified sample. An example of labeled DNA run on an agarose gel is shown in Figure 3.

**Step 98**

A sample region of a scanned array is shown in Figure 4. Arrays should have a clean and uniform background, with little particulate debris or dust. Features should be evenly arranged and have a uniform shape. Typically, we label immunoprecipitated material with Cy5 and input DNA material with Cy3, and set array visualization so that Cy5 signals appear red and Cy3 signals appear green. Thus, for a successful immunoprecipitation, we expect a subpopulation of red spots, indicating features that have been enriched in the immunoprecipitation channel. We expect few, if any, non-control features that appear bright green, indicating features enriched in the input DNA channel. We expect the majority of spots to be yellow or yellowish-green, indicating no enrichment. Note that this assumes there are relatively few binding events and that the array has been scanned to balance the distributions of signal intensities in the two channels.

**Step 104**

Figure 5 shows a progression of data from array signal to identification of genomic regions that have been enriched for binding to verification by gene-specific PCR. The array images (Fig. 5a,d) are from slides that were part of a multislide, whole-genome array that detects binding of the initiating form of RNA polymerase II. The slides were scanned and processed, as described, to generate ratios of signal intensities. The subsequent ratios were processed to identify significantly enriched features and binding events.

When we plot enrichment ratios against chromosomal position of array features, we expect binding events to have two basic characteristics (Fig. 5b). Features identifying a binding event should show significantly higher enrichment ratios than background. Observed enrichment ratios generally range from 2-fold to 15-fold, although ratios will vary based on antibody performance and can be as high as 80-fold in some cases. Binding events should be supported by multiple features: two or more adjacent features should both indicate enrichment above background. The number of probes that are expected to identify a binding event will depend on the degree of sonication, which influences the size of DNA fragments in the ChIP and the tiling density of the array, which affects the average number of probes that could reasonably detect immuno-enriched fragments, given the distribution of fragment sizes and whether the protein binds over large or small regions.

Legitimate binding events can also be confirmed by PCR (Fig. 5c). PCR should be performed using primers that are designed to amplify the genomic regions that are bound and run the products should be run on an agarose gel. The signals that are generated using the immunoprecipitated DNA sample should be compared to the signal generated using the input DNA sample. There should be more signal in the immunoprecipitated DNA sample when equivalent amounts of immunoprecipitated DNA and input DNA are used in the PCR. A serial dilution of the input DNA PCR reactions should be performed to ensure that all of the reactions are in the linear range of detection.

**Acknowledgments**

The protocol and suggestions presented here represent the experience accumulated over several years of experimentation by a number of researchers in the lab. In particular, we would like to acknowledge Dmitry




Pokholok and Brett Chevalier for their contributions to streamlining the protocol. Other particularly useful contributions or observations were offered by Duncan Odom, Beth Jacobsen, Richard Jenner, Laurie Boyer, Matthew Guenther, Roshan Kumar, Heather Plasterer, David Reynolds and Jane Yoo. We thank current members of the lab, particularly Megan Cole, Matthew Guenther, Beth Jacobsen, Richard Jenner and Michael Kagey, for critical reading of this manuscript.

## References

1. Solomon MJ, Larsen PL, Varshavsky A. Mapping protein-DNA interactions *in vivo* with formaldehyde: evidence that histone H4 is retained on a highly transcribed gene. *Cell* 1988;53:937–947. [PubMed: 2454748]
2. Dedon PC, Soultis JA, Allis CD, Gorovsky MA. Formaldehyde cross-linking and immunoprecipitation demonstrate developmental changes in H1 association with transcriptionally active genes. *Mol. Cell. Biol* 1991;11:1729–1733. [PubMed: 1996118]
3. Orlando V, Paro R. Mapping Polycomb-repressed domains in the bithorax complex using *in vivo* formaldehyde cross-linked chromatin. *Cell* 1993;75:1187–1198. [PubMed: 7903220]
4. Reid JL, Iyer VR, Brown PO, Struhl K. Coordinate regulation of yeast ribosomal protein genes is associated with targeted recruitment of Esa1 histone acetylase. *Mol. Cell* 2000;6:1297–1307. [PubMed: 11163204]
5. Ren B, et al. Genome-wide location and function of DNA binding proteins. *Science* 2000;290:2306–2309. [PubMed: 11125145]
6. Iyer VR, et al. Genomic binding sites of the yeast cell-cycle transcription factors SBF and MBF. *Nature* 2001;409:533–538. [PubMed: 11206552]
7. Lieb JD, Liu X, Botstein D, Brown PO. Promoter-specific binding of Rap1 revealed by genome-wide maps of protein-DNA association. *Nature Genet* 2001;28:327–334. [PubMed: 11455386]
8. Kurdستاني SK, Robyr D, Tavazoie S, Grunstein M. Genome-wide binding map of the histone deacetylase Rpd3 in yeast. *Nature Genet* 2002;31:248–254. [PubMed: 12089521]
9. Lee TI, et al. Transcriptional regulatory networks in *Saccharomyces cerevisiae*. *Science* 2002;298:799–804. [PubMed: 12399584]
10. Ren B, et al. E2F integrates cell cycle progression with DNA repair, replication, and G(2)/M checkpoints. *Genes Dev* 2002;16:245–256. [PubMed: 11799067]
11. Weinmann AS, Yan PS, Oberley MJ, Huang TH, Farnham PJ. Isolating human transcription factor targets by coupling chromatin immunoprecipitation and CpG island microarray analysis. *Genes Dev* 2002;16:235–244. [PubMed: 11799066]
12. Kim TH, et al. A high-resolution map of active promoters in the human genome. *Nature* 2005;436:876–880. [PubMed: 15988478]
13. Pokholok DK, et al. Genome-wide map of nucleosome acetylation and methylation in yeast. *Cell* 2005;122:517–527. [PubMed: 16122420]
14. Lee TI, et al. Control of developmental regulators by Polycomb in human embryonic stem cells. *Cell* 2006;125:301–313. [PubMed: 16630818]
15. Nal B, Mohr E, Ferrier P. Location analysis of DNA-bound proteins at the whole-genome level: untangling transcriptional regulatory networks. *Bioessays* 2001;23:473–476. [PubMed: 11385626]
16. Hanlon SE, Lieb JD. Progress and challenges in profiling the dynamics of chromatin and transcription factor binding with DNA microarrays. *Curr. Opin. Genet. Dev* 2004;14:697–705. [PubMed: 15531167]
17. Blais A, Dynlacht BD. Constructing transcriptional regulatory networks. *Genes Dev* 2005;19:1499–1511. [PubMed: 15998805]
18. Sikder D, Kodadek T. Genomic studies of transcription factor-DNA interactions. *Curr. Opin. Chem. Biol* 2005;9:38–45. [PubMed: 15701451]
19. Carter NP, Vetrie D. Applications of genomic microarrays to explore human chromosome structure and function. *Hum. Mol. Genet* 2004;13(Spec No 2):R297–R302. [PubMed: 15358737]
20. van Steensel B, Henikoff S. Epigenomic profiling using microarrays. *Biotechniques* 2003;35:346–357. [PubMed: 12951776]


21. MacAlpine DM, Bell SP. A genomic view of eukaryotic DNA replication. *Chromosome Res* 2005;13:309–326. [PubMed: 15868424]
22. Woodfine K, Carter NP, Dunham I, Fiegler H. Investigating chromosome organization with genomic microarrays. *Chromosome Res* 2005;13:249–257. [PubMed: 15868419]
23. Lei EP, Krebber H, Silver PA. Messenger RNAs are recruited for nuclear export during transcription. *Genes Dev* 2001;15:1771–1782. [PubMed: 11459827]
24. Ishii K, Arib G, Lin C, Van Houwe G, Laemmli UK. Chromatin boundaries in budding yeast: the nuclear pore connection. *Cell* 2002;109:551–562. [PubMed: 12062099]
25. Rodriguez-Navarro S, et al. Sus1, a functional component of the SAGA histone acetylase complex and the nuclear pore-associated mRNA export machinery. *Cell* 2004;116:75–86. [PubMed: 14718168]
26. Takizawa PA, DeRisi JL, Wilhelm JE, Vale RD. Plasma membrane compartmentalization in yeast by messenger RNA transport and a septin diffusion barrier. *Science* 2000;290:341–344. [PubMed: 11030653]
27. Tenenbaum SA, Carson CC, Lager PJ, Keene JD. Identifying mRNA subsets in messenger ribonucleoprotein complexes by using cDNA arrays. *Proc. Natl. Acad. Sci. USA* 2000;97:14085–14090. [PubMed: 11121017]
28. Brown V, et al. Microarray identification of FMRP-associated brain mRNAs and altered mRNA translational profiles in fragile X syndrome. *Cell* 2001;107:477–487. [PubMed: 11719188]
29. Hieronymus H, Silver PA. Genome-wide analysis of RNA-protein interactions illustrates specificity of the mRNA export machinery. *Nature Genet* 2003;33:155–161. [PubMed: 12524544]
30. Shepard KA, et al. Widespread cytoplasmic mRNA transport in yeast: identification of 22 bud-localized transcripts using DNA microarray analysis. *Proc. Natl Acad. Sci. USA* 2003;100:11429–11434. [PubMed: 13679573]
31. Sanchez-Elsner T, Gou D, Kremmer E, Sauer F. Noncoding RNAs of trithorax response elements recruit *Drosophila* Ash1 to Ultrabithorax. *Science* 2006;311:1118–1123. [PubMed: 16497925]
32. Orlando V, Strutt H, Paro R. Analysis of chromatin structure by *in vivo* formaldehyde cross-linking. *Methods* 1997;11:205–214. [PubMed: 8993033]
33. Kuo MH, Allis CD. *In vivo* cross-linking and immunoprecipitation for studying dynamic Protein:DNA associations in a chromatin environment. *Methods* 1999;19:425–433. [PubMed: 10579938]
34. Buck MJ, Lieb JD. ChIP-chip: considerations for the design, analysis, and application of genome-wide chromatin immunoprecipitation experiments. *Genomics* 2004;83:349–360. [PubMed: 14986705]
35. Boyer LA, et al. Core transcriptional regulatory circuitry in human embryonic stem cells. *Cell* 2005;122:947–956. [PubMed: 16153702]
36. Boyer LA, et al. Polycomb complexes repress developmental regulators in murine embryonic stem cells. *Nature* 2006;441:349–353. [PubMed: 16625203]
37. Lieb JD. Genome-wide mapping of protein-DNA interactions by chromatin immunoprecipitation and DNA microarray hybridization. *Methods Mol. Biol* 2003;224:99–109. [PubMed: 12710669]
38. Oberley MJ, Farnham PJ. Probing chromatin immunoprecipitates with CpG-island microarrays to identify genomic sites occupied by DNA-binding proteins. *Methods Enzymol* 2003;371:577–596. [PubMed: 14712730]
39. Bannister AJ, Kouzarides T. Histone methylation: recognizing the methyl mark. *Methods Enzymol* 2004;376:269–288. [PubMed: 14975312]
40. Bernstein BE, Humphrey EL, Liu CL, Schreiber SL. The use of chromatin immunoprecipitation assays in genome-wide analyses of histone modifications. *Methods Enzymol* 2004;376:349–360. [PubMed: 14975317]
41. Chaya D, Zaret KS. Sequential chromatin immunoprecipitation from animal tissues. *Methods Enzymol* 2004;376:361–372. [PubMed: 14975318]
42. Ciccone DN, Morshead KB, Oettinger MA. Chromatin immunoprecipitation in the analysis of large chromatin domains across murine antigen receptor loci. *Methods Enzymol* 2004;376:334–348. [PubMed: 14975316]

43. Oberley MJ, Tsao J, Yau P, Farnham PJ. High-throughput screening of chromatin immunoprecipitates using CpG-island microarrays. *Methods Enzymol* 2004;376:315–334. [PubMed: 14975315]
44. Ren B, Dynlacht BD. Use of chromatin immunoprecipitation assays in genome-wide location analysis of mammalian transcription factors. *Methods Enzymol* 2004;376:304–315. [PubMed: 14975314]
45. Bernstein BE, Liu CL, Humphrey EL, Perlstein EO, Schreiber SL. Global nucleosome occupancy in yeast. *Genome Biol* 2004;5:R62. [PubMed: 15345046]
46. Cawley S, et al. Unbiased mapping of transcription factor binding sites along human chromosomes 21 and 22 points to widespread regulation of noncoding RNAs. *Cell* 2004;116:499–509. [PubMed: 14980218]
47. Liu CL, et al. Single-nucleosome mapping of histone modifications in *S. cerevisiae*. *PLoS Biol* 2005;3:e328. [PubMed: 16122352]
48. Rada-Iglesias A, et al. Binding sites for metabolic disease related transcription factors inferred at base pair resolution by chromatin immunoprecipitation and genomic microarrays. *Hum. Mol. Genet* 2005;14:3435–3447. [PubMed: 16221759]
49. Rao B, Shibata Y, Strahl BD, Lieb JD. Dimethylation of histone H3 at lysine 36 demarcates regulatory and nonregulatory chromatin genome-wide. *Mol. Cell. Biol* 2005;25:9447–9459. [PubMed: 16227595]
50. Yuan GC, et al. Genome-scale identification of nucleosome positions in *S. cerevisiae*. *Science* 2005;309:626–630. [PubMed: 15961632]
51. Bieda M, Xu X, Singer MA, Green R, Farnham PJ. Unbiased location analysis of E2F1-binding sites suggests a widespread role for E2F1 in the human genome. *Genome Res* 2006;16:595–605. [PubMed: 16606705]
52. Hughes TR, et al. Functional discovery via a compendium of expression profiles. *Cell* 2000;102:109–126. [PubMed: 10929718]

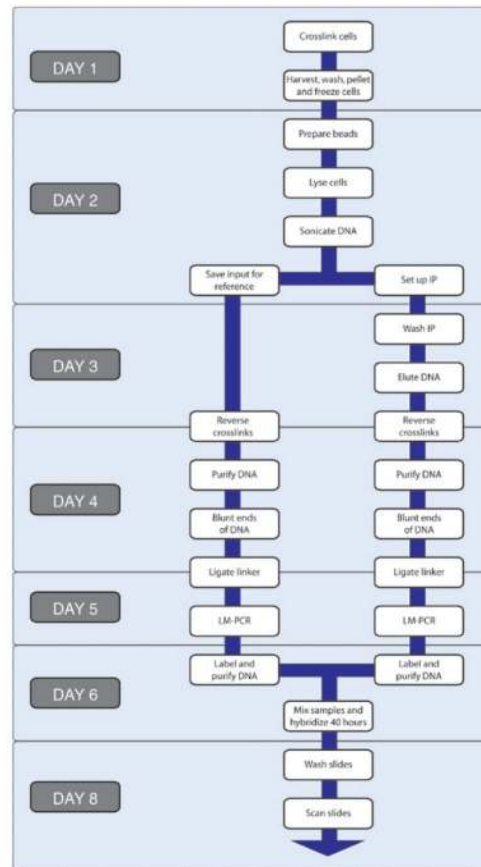


**Figure 1.**

A sample timeline for the ChIP-Chip protocol. Individual steps are shown in white boxes. Steps that are typically performed on the same day are grouped by day, which is indicated in gray boxes. IP, immunoprecipitation; LM-PCR, ligation-mediated PCR.

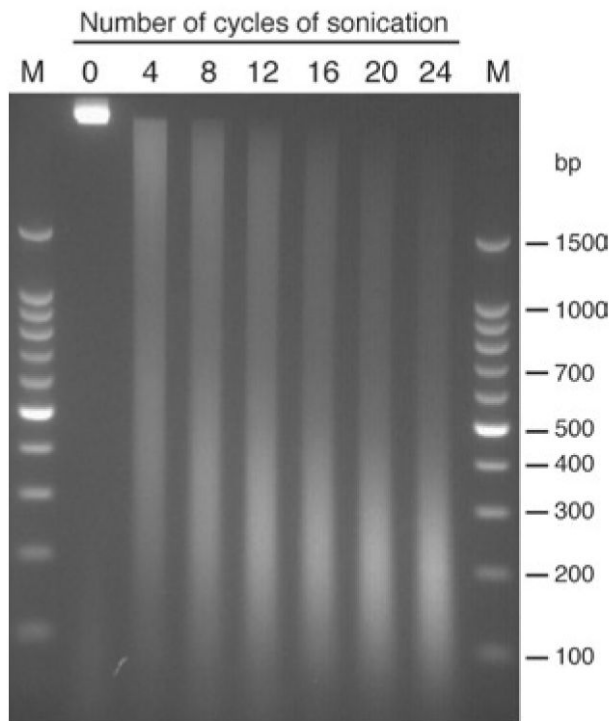


**Figure 2.**

Results of varying degrees of sonication on fragment size. A total of  $2 \times 10^8$  crosslinked Jurkat cells were sonicated using the following conditions: Misonix 3000 sonicator with microtip; power 7; 24 cycles (30 s sonication, 90 s rest). Samples were removed at various times, crosslinks were reversed and DNA-purified, and run on a 2% agarose gel. Lanes with molecular weights are labeled M and lanes with sonicated material are labeled with the number of cycles of sonication. Sizes of molecular-mass markers are indicated. Twelve cycles of sonication provide a good degree of sonication. Four cycles of sonication results in undersonicated DNA. Note that sonication beyond 16 cycles results in little change in fragment size and is likely to result in oversonicated material.

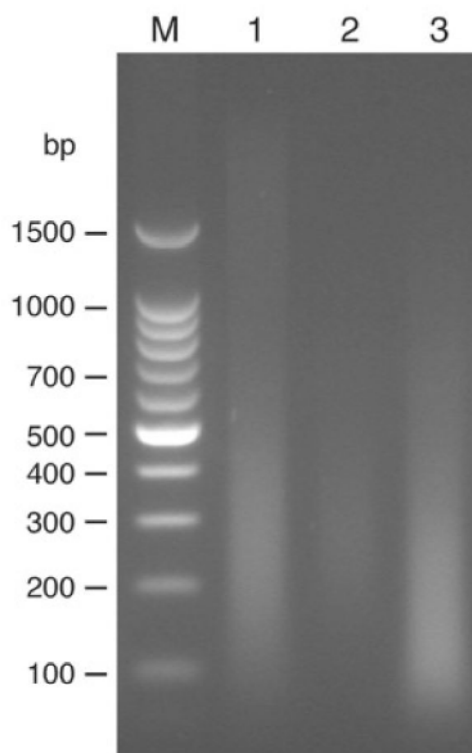


**Figure 3.** 2% agarose gel showing an example of input DNA, DNA after LMPCR amplification and labeled DNA. The lane with molecular-mass markers is labeled M. Input DNA is in Lane 1, LMPCR DNA is in Lane 2, labeled DNA is in Lane 3. Approximately 200–300 ng of DNA are loaded. Note that labeled DNA will fluoresce (and thus appear more abundant) when photographed with UV illumination due to the presence of dye.

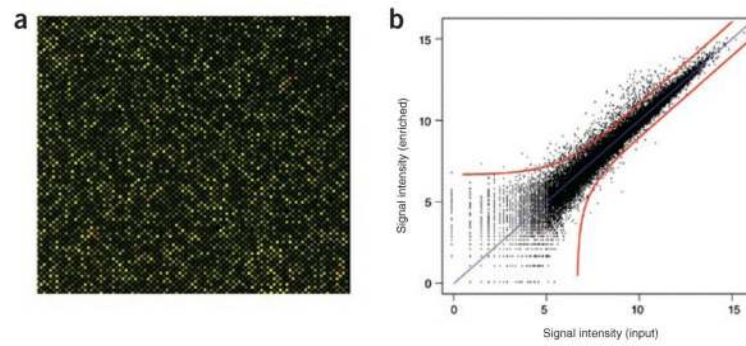

**Figure 4.**

Samples of hybridized arrays and scatterplots. **(a)** A portion of a hybridized and scanned array. The material is enriched for RNA polymerase II binding. Red features indicate enrichment. Bright green features are hybridization controls. Yellow and yellow-green features indicate no enrichment. **(b)** A plot showing processed results of the hybridization and scan shown in **a**. Signal intensities (log base 2) for the RNA polymerase II-enriched DNA are shown on the *y*-axis. Signal intensities (log base 2) for the input DNA are shown on the *x*-axis. Features outside the boundaries, indicated by the red lines, are significantly enriched at a threshold of  $P < 0.001$ . Features that are below the diagonal and appear significantly enriched are usually control features or dust particles that have been mistakenly identified as features. Features that are above the diagonal and appear significantly enriched usually represent immuno-enriched features and indicate protein binding.



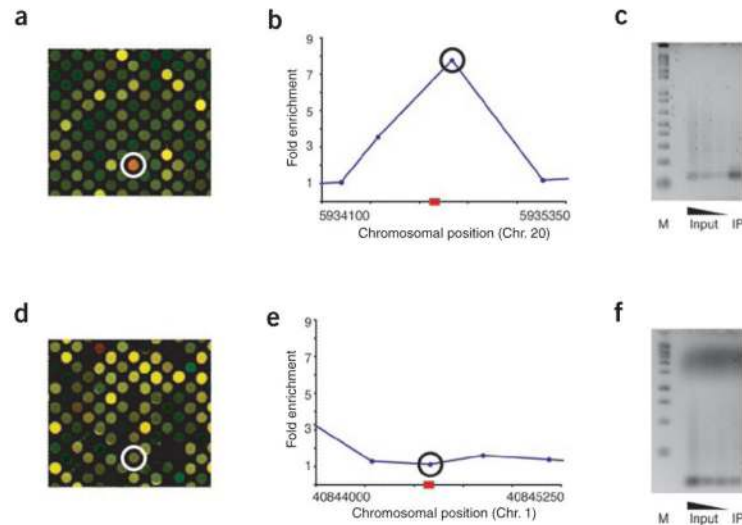

**Figure 5.**

Examples of data from an array that are processed to identify genomic regions that are enriched for binding and verification by gene-specific PCR. **(a)** Another portion of the hybridized and scanned array shown in Figure 4. The material is enriched for RNA polymerase II binding. The circled feature is red, indicating enrichment. Yellow and yellow-green features indicate no enrichment. **(b)** A plot showing ratios of the RNA polymerase II enrichment. Enrichment is plotted on the y-axis. Genomic location is plotted on the x-axis. Probe location is shown by the blue circles. The circled probe represents the ratio seen for the circled feature shown in **a**. The red box indicates the region checked using gene-specific PCR to confirm enrichment in **c**. **(c)** Gene-specific PCR, confirming enrichment. Primers were used to PCR the genomic region, indicated by the red box in **b** for immunoprecipitated DNA and dilutions of input DNA. Molecular-mass markers are shown in the left-most lane, which is labeled M. Three-fold dilutions of input DNA are shown in lanes 1–3, as indicated. The immunoprecipitation (IP) is shown in the last lane. The amount of immunoprecipitate that is added to the PCR is equivalent to the last dilution of input DNA. **(d)** A portion of a hybridized and scanned array. The hybridized material is enriched for RNA polymerase II binding. The circled feature shows no enrichment. Yellow and yellow-green features indicate no enrichment. **(e)** A plot showing ratios of the RNA polymerase II enrichment. Enrichment is plotted on the y-axis. Genomic location is plotted on the x-axis. Probe location is shown by the blue circles. The circled probe represents the ratio seen for the circled feature shown in **d**. The red box indicates the region checked using gene-specific PCR to confirm enrichment in **f**. **(f)** Gene-specific PCR confirming lack of enrichment. Primers were used to PCR the genomic region, indicated by the red box in **d** for immunoprecipitated DNA and dilutions of input DNA. Molecular-mass markers are shown in the left-most lane, which is labeled M. Three-fold dilutions of input DNA are shown in lanes 1–3, as indicated. The immunoprecipitation is shown in the last lane. The amount of immunoprecipitate added to the PCR is equivalent to the last dilution of input DNA.

TABLE 1

Troubleshooting table.

PROBLEM	POSSIBLE REASON	SOLUTION
Insufficiently sheared chromatin (detected at Step 42)	Over-crosslinked material (Step 1)	Reduce crosslinking (Step 1)
	Solution foamed (Step 17)	Keep sonication tip centered and within 1 cm of bottom of tube (Step 15)
	Solution foamed (Step 17)	If the solution has foamed before, make sure that all bubbles are removed (Step 17)
Low signal:noise or failure to immunoprecipitated DNA (detected at Step 42 or Step 103)	Under-sonicated (Step 17)	Increase sonication (either power settings or number of rounds of sonication) (Step 17)
	Over-crosslinked material (epitope obscured or denatured) (Step 1)	Reduce crosslinking time (Step 1)
	Under-crosslinked (protein insufficiently crosslinked to DNA) (Step 1)	Increase crosslinking time. We have gone as long as 16 h at 4 °C (Step 1)
	Ineffective antibody (Step 7)	If possible, try a different antibody (Step 7)
Debris on slides (detected at Step 98)	Poor conditions for immunoprecipitation (Step 21)	Change salt concentrations of immunoprecipitation and subsequent washes (Step 21 and Step 25)
	Particulates not removed from labeled material (Step 76, Step 86)	Spin labeled DNA at high speed to remove particulates (Step 76, Step 86)
	Particulates during wash (Step 91, Step 92, Step 93, Step 94, Step 95)	Use clean glassware and filter necessary wash solutions. If using an ozone-scavenger solution, dissolve solute completely (Step 89, Step 90)
Noise in microarray analysis (abnormally wide distribution of immunoprecipitated/input DNA signal ratios) (detected at Step 103)	Lines/streaks on slides (Step 98)	Transfer slides quickly between washes to minimize drying. When removing slides from final wash, avoid sudden movements (Step 93, Step 95, Step 97)
	Over-sonicated material (Step 17)	Reduce sonication (either power settings or number of rounds of sonication) (Step 17)
	Beads overwashed (Step 25)	Minimize wash times; reduce number of washes (Step 25)
	Incubation period too long for crosslink reversal (Step 31, Step 32)	Reduce incubation for crosslink reversal to 6–8 h (Step 31, Step 32)
	DNA degraded by phenol (Step 36, Step 47)	Replace P:C:IA solution (Step 36, Step 47)
	Poor recovery of DNA (Step 40, Step 50, Step 56, Step 63)	Take care when resuspending pellets (Step 40, Step 50, Step 56, Step 63)
	Poor blunting (Step 45)	Keep reaction temperature low and incubation time short (Step 45)
Small DNA fragments post-labeling (Step 77)	Use a clean-up method that limits purification of small (60 bp or less) fragments (Step 75)	