**Research Paper** 

# Clinical mutational profiling of 1006 lung cancers by next generation sequencing

Peter B. Illei<sup>1,\*</sup>, Deborah Belchis<sup>1,\*</sup>, Li-Hui Tseng<sup>1,2</sup>, Doreen Nguyen<sup>1</sup>, Federico De Marchi<sup>1,3</sup>, Lisa Haley<sup>1</sup>, Stacy Riel<sup>1</sup>, Katie Beierl<sup>1</sup>, Gang Zheng<sup>1</sup>, Julie R. Brahmer<sup>4</sup>, Frederic B. Askin<sup>1</sup>, Christopher D. Gocke<sup>1,4</sup>, James R. Eshleman<sup>1,4</sup>, Patrick M. Forde<sup>4</sup> and Ming-Tseh Lin<sup>1</sup>

<sup>1</sup>Department of Pathology, Johns Hopkins University School of Medicine, Johns Hopkins Hospital, Baltimore, Maryland, USA

<sup>2</sup>Department of Medical Genetics, National Taiwan University Hospital, Taipei, Taiwan

<sup>3</sup>Division of Hematology and Bone Marrow Transplantation, University of Udine Hospital, Udine, Italy

<sup>4</sup>Department of Oncology, Johns Hopkins University School of Medicine, Johns Hopkins Hospital, Baltimore, Maryland, USA

\*These authors contributed equally to this work

Correspondence to: Peter B. Illei, email: mlin36@jhmi.edu Ming-Tseh Lin, email: mlin36@jhmi.edu

Keywords: lung, mutation, sequencing, cancer, profiling

Received: February 21, 2017 Accepted: May 10, 2017 Published: May 20, 2017

Copyright: Illei et al. This is an open-access article distributed under the terms of the Creative Commons Attribution License 3.0 (CC BY 3.0), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited.

#### ABSTRACT

Analysis of lung adenocarcinomas for actionable mutations has become standard of care. Here, we report our experience using next generation sequencing (NGS) to examine AKT1, BRAF, EGFR, ERBB2, KRAS, NRAS, and PIK3CA genes in 1006 non-small cell lung cancers in a clinical diagnostic setting. NGS demonstrated high sensitivity. Among 760 mutations detected, the variant allele frequency (VAF) was 2–5% in 33 (4.3%) mutations and 2–10% in 101 (13%) mutations. A single bioinformatics pipeline using Torrent Variant Caller, however, missed a variety of EGFR mutations. Mutations were detected in KRAS (36% of tumors), EGFR (19%) including 8 (0.8%) within the extracellular domain (4 at codons 108 and 4 at codon 289), BRAF (6.3%), and PIK3CA (3.7%). With a broader reportable range, exon 19 deletion and p.L858R accounted for only 36% and 26% of EGFR mutations and p.V600E accounted for only 24% of **BRAF** mutations. NGS provided accurate sequencing of complex mutations seen in 19% of EGFR exon 19 deletion mutations. Doublet (compound) EGFR mutations were observed in 29 (16%) of 187 EGFR-mutated tumors, including 69% with two non-p. L858R missense mutations and 24% with p.L858 and non-p.L858R missense mutations. Concordant VAFs suggests doublet EGFR mutations were present in a dominant clone and cooperated in oncogenesis. Mutants with predicted impaired kinase, observed in 25% of BRAF-mutated tumors, were associated with a higher incidence of concomitant activating KRAS mutations. NGS demonstrates high analytic sensitivity, broad reportable range, quantitative VAF measurement, single molecule sequencing to resolve complex deletion mutations, and simultaneous detection of concomitant mutations.

#### **INTRODUCTION**

Approximately 30–40% of Asian patients and 10–15% of Caucasian patients with lung adenocarcinoma harbor activating mutations in the *epidermal growth factor* 

*receptor* (*EGFR*) gene. Gefitinib, erlotinib and afatinib are tyrosine kinase inhibitors (TKI) approved by the Food and Drug Administration (FDA) of the United States for treatment of patients with *EGFR*-mutated lung cancers [1–3]. In 2011, a provisional clinical opinion from the American Society of Clinical Oncology recommended testing for *EGFR* mutations in patients with metastatic lung cancer to predict response to TKI therapy [4]. Molecular testing guidelines for selection of lung cancer patients for TKI therapy have been published and are currently under revision by the College of American Pathologists, International Association for the Study of Lung Cancer, and Association for Molecular Pathology [5].

A variety of molecular diagnostic assays have been clinically validated for detection of EGFR mutations [6]. Although the prior gold standard of Sanger sequencing covers all EGFR mutations within exons 18-21, it's analytic sensitivity (20-40% tumor cellularity) may not be adequate in the clinical diagnostic setting where specimens containing low tumor cellularity are not uncommon [7, 8]. The analytic sensitivity can be improved to approximately 5% variant allele frequency (VAF) (10% tumor cellularity) using pyrosequencing, 1% VAF using mutation-specific real time PCR assays, or even less than 1% by droplet digital PCR. Currently, there are two assays approved by the FDA for testing EGFR mutations in lung cancers, the cobas EGFR mutation test (Roche Molecular Systems, Branchburg, NJ) and the therascreen EGFR RGQ PCR Kit (Qiagen, Hilden, Germany) [9-13]. Both assays detect hot spot EGFR mutations by multiple separate runs of mutation-specific real-time PCR assays. These assays are not able to detect less common mutations outside the reportable ranges. A total of 150 ng DNA is needed for the cobas test. DNA input has not been quantified for the therascreen test. We have shown that 44% of specimens submitted for clinical mutational profiling were taken by biopsy or fine needle aspiration [14]. DNA isolated from biopsy or fine needle aspiration specimens containing limited tissue may not be sufficient.

Multiplexed genotyping platforms are replacing the traditional "one test-one drug" paradigm not only because of continuous expansion of predictive markers for targeted therapeutics but also often limited tissues submitted to the clinical diagnostic laboratories. Primer extension-based assays with a multiplex design, such as the Sequenom MassARRAY system, detect multiple hotspots within a panel of genes including EGFR in a single reaction while retaining an analytic sensitivity of 5% or less VAF [15, 16]. In the era of precision cancer medicine, molecular diagnostics assays with a higher analytic sensitivity and a broader reportable range are warranted to provide a more comprehensive mutational profiling. Next generation sequencing (NGS) technology has led to a revolution in genome discovery and will soon become the most costeffective multiplexed sequencing platform in the setting of clinical care [7, 17]. We have previously validated a NGS platform using the AmpliSeq Cancer Hotspot Panel and Personal Genome Machine in a Clinical Laboratory Improvement Amendments (CLIA)-certified laboratory [18]. In this retrospective analysis for quality assessment, we surveyed our experience with clinical mutation detection of *AKT1*, *BRAF*, *ERBB2*, *EGFR*, *KRAS*, *NRAS*, and *PIK3CA* genes in 1006 lung cancer specimens using this NGS assay, including false negative calling, the capability of detecting complex deletion mutations within exon 19 of the *EGFR* gene, a high frequency of doublet (compound) *EGFR* mutations with concordant VAFs, and the association of kinase impaired *BRAF* mutations with activating *RAS* mutations.

# RESULTS

# Positive and negative controls

No mutations were detected in 88 runs of the negative control specimen while all mutations in the positive control specimens were detected. The observed VAFs in the positive controls were highly consistent throughout the test period, demonstrating that NGS is a precise quantification assay for VAF (Supplementary Table 1).

# Mutations missed by Torrent Variant Caller

Specimens with prior TKI therapy were not included in this study. One mutation was detected in 560 of 1006 tumors, 2 mutations in 70 tumors, and 3 mutations in 3 tumors. Our analysis pipeline included both Torrent Variant Caller and direct visual inspection of all amplicons within the reportable range using Integrative Genomics Viewer (IGV). A total of 15 mutations detected by IGV inspection were missed by Torrent Variant Caller (Table 1). These included 12 EGFR mutations (3 missense mutations of exon 18, 3 deletion mutations of exon 19, 3 insertion/duplication mutations of exon 20 and 3 p.L858R mutations of exon 21), 2 PIK3CA missense mutations and one ERBB2 duplication mutation. VAFs ranged from 3.2% to 65% (Figure 1A). Eleven of the 15 false negative calls by Torrent Variant Caller occurred within the first year. The false negative calls were most likely related to bioinformatics pipelines of the Torrent Variant Caller, which has been improved with the updated versions. The recent Torrent Variant Caller did not miss mutations in lung cancers. All insertion/deletion (indel) mutations detected within EGFR exons 19 and 20 by a capillary electrophoresis sizing assay were also detected by IGV inspection. Accordingly, sizing by capillary electrophoresis was discontinued in April 2014.

# EGFR mutations

A total of 49 unique *EGFR* variants/mutations were detected in 187 (19%) tumors, including 29 tumors with 2 *EGFR* mutations and 10 variants/mutations outside exons 18–21 involving extracellular domain of EGFR (four codon 108 mutations, four codon 289 mutations and two possible germline variants, p.E282K and p.H584R) (Table 2 and Supplementary Table 2). All except p.L858fs\*45 caused

Gene <sup>a</sup>	cDNA change	Amino acid change	VAF <sup>b</sup>	VC version
EGFR (Re)	c.2573T>G	p.L858R	8.1	3.4.51874
EGFR (Re)	c.2171G>C	p.G724A	16	3.4.51874
PIK3CA (Bx)	c.1634A>C	p.E545A	10	3.4.51874
EGFR (FNA)	c.2237_2255delinsT <sup>d</sup>	p.E746_S752delinsV	3.2	3.4.51874
EGFR (Bx)	c.2156G>C	p.G719A	8.9	3.4.51874
EGFR (Re)	c.2573T>G	p.L858R	6.1	3.4.51874
EGFR (Re)	c.2573T>G	p.L858R	12	3.4.51874
EGFR (FNA)	c.2236_2250del	p.E746_A750del	5.5	3.6.63335
EGFR (Bx)	c.2300_2308dup	p.A767_V769dup	9.4	3.6.63335
EGFR (Bx)	c.2302_2303insCGCTGGCCA	p.A767_S768insTLA	5	3.6.63335
EGFR (Re)	c.2156G>C	p.G719A	32	3.6.63335
ERBB2 (FNA)	c.2313_2324dup	p.A771_M774dup	7.6	3.6.63335
EGFR (Bx)	c.2311_2312delinsGGGTT <sup>e</sup>	p.N771delinsGF	65	3.6.63335
PIK3CA (Bx)	c.1193G>T	p.R398L	5.5	3.6.63335
EGFR (Re)	c.2239_2256del	p.L747_S752del	36	3.6.63335

Table 1: Mutations detected by integrative genomics viewer inspection but missed by Torrent Variant Caller (VC)

<sup>a</sup>Re: resection specimens, Bx: biopsy specimens; FNA: fine needle aspiration specimens. Two specimens (*EGFR* p.N771delinsGF mutation and *PIK3CA* p.R398L mutation) were concentrated by using the Amicon filter.

<sup>b</sup>VAF: variant allelic frequency. Depth of coverage was more than 1000 reads except the specimen with *PIK3CA* p.R398 mutation (450 reads).

<sup>c</sup>The current VC version (5.0.2.1) did not miss mutations detected by IGV since December 2015.

<sup>d</sup>The c.2237\_2255delinsT mutation was composed of c.2237\_2254del and c.2255C>T. Both were missed by Torrent Variant Caller.

<sup>e</sup>The c.2311\_2312delinsGGGTT mutation was composed of c.2310-2311insGGG, c.2311A>T and c.2312A>T. The Torrent Variant Caller detected c.2311A>T and c.2312A>T, but not c.2310-2311insGGG.

	Lung cancer		Variant allele freque	ency <sup>a</sup>
	(n = 1006)	≤5%	$\leq 10\%$	≤ 20%
AKT1	4 (0.4%)	NA	NA	NA
BRAF	63 (6.3%)	7.7%	25% ( <i>P</i> = 0.009)	63% ( <i>P</i> < 0.001)
EGFR	187 (19%)	2.2%	7.8% (P = 0.07)	27% ( <i>P</i> = 0.003)
ERBB2	13 (1.3%)	NA	NA	NA
KRAS	362 (36%)	4.1%	13%	39%
NRAS	11 (1.1%)	NA	NA	NA
PIK3CA	37 (3.7%)	13%	37% ( <i>P</i> < 0.001)	63% ( <i>P</i> = 0.003)

Table 2: Mutational profiling of 1006 lung cancers and variant allele frequency

NA: not analyzed.

<sup>a</sup>Percentage of mutations with variant allele frequencies of  $\leq 5\%$ ,  $\leq 10\%$  and  $\leq 20\%$ . Data of all 1074 specimens were included for analysis of variant allele frequencies. *P* values were analyzed by comparing with *KRAS* variant allele frequencies.

either missense mutations or in-frame indel mutations. Distribution of the 216 *EGFR* mutations was shown in Supplementary Figure 1 and frequencies of TKI sensitive or resistant mutations among the 187 *EGFR*-mutated tumors were shown in Supplementary Figure 2.

At the level of cDNA, there were 5 unique simple exon 19 deletion mutations (defined as an exon 19 mutation with no additional nucleotide changes within 1–2 adjacent codons) and 11 unique complex exon 19 deletion mutations (defined as an exon 19 deletion accompanied with one or more nucleotide changes in *cis* position within 1–2 adjacent codons) (Table 3). Complex mutations were seen in 15 (19%) of 78 exon 19 deletion mutations, including 12 with one deletion and 1 or 2 single nucleotide changes, 2 with two deletions, and one with two deletions and one single nucleotide change (Figure 1B). All the accompanied single nucleotide changes were located 3' to the deletion. All 5 simple deletion mutations were observed in 2 or more tumors while most complex deletion mutations were observed in one tumor.

# Mutations in the mitogen-activated protein kinase (MAPK) pathway

A total of 19 unique *KRAS* mutations were detected in 362 (36%) tumors, including 3 tumors with two *KRAS* mutations (Table 2 and Supplementary Table 3). These mutations were located within exon 2 in 326 (89%) of 365 *KRAS* mutations, exon 3 in 9.0% mutations, and exon 4 in 1.6% mutations. The two most common *KRAS* mutations were p.G12C (34%) and p.G12V (22%). Mutations located outside codon 12 or 13 were observed in 43 tumors (12% of *KRAS*-mutated tumors). Six unique *NRAS* mutations were detected in 11 tumors (1.1%), 6 at codon 12, 4 at codon 61 and one p.G48A (Table 2 and Supplementary Table 3).

Twenty-four unique *BRAF* mutations were detected in 63 (6.3%) tumors (Table 2 and Supplementary Table 4). The most common residue involved by *BRAF* mutations

A (3.2%)

was codon 600 (25%), followed by codon 469 (18%), codon 466 (16%), codon 594 (14%), and codon 601 (11%) (Supplementary Figure 3). The p.V600E mutant accounted for only 24% of *BRAF* mutations and 40% of *BRAF* mutations occurred within exon 11. Of the 63 *BRAF*-mutated tumors, BRAF kinase activity was predicted to be elevated in 33 (52%), impaired in 16 (25%) and unknown in 14 (22%) tumors. The most common kinase-impaired *BRAF* mutants are p.G466V (5) and p.D594G (5), followed by p.D594N (3), p.G466E (1), p.G466R (1) and p.D594H (1).

# Mutations in the phosphatidylinositol 3-kinase (mTOR) pathway


*AKT1* mutations were detected in 3 tumors with p.E17K and one with p.R25H (Table 2 and Supplementary Table 5). Seventeen unique *PIK3CA* mutations were detected in 37 (3.7%) tumors (Table 2 and Supplementary Table 5). All except p.R401\* were missense mutations. Most mutations were located within exons 9 (57%), 20 (22%) and 4 (11%). The 3 most common codons (E542, E545 and H1047) accounted for only 65% of the *PIK3CA* mutations.

# ERBB2 mutations

Five unique *ERBB2* mutations were detected in 13 (1.3%) tumors (Table 2 and Supplemental Table 5). There were 12 insertion/duplication mutations within exon 20 and one missense mutation within exon 19. The most common mutation was p.A771\_M774dup (62% of *ERBB2* mutations).

# Variant allele frequency

A total of 760 mutations were detected in the 1074 specimens. VAFs were 2-5% in 33 (4.3%) mutations, 2-10% in 101 (13%) mutations, and 2-20% in 289 (38%)

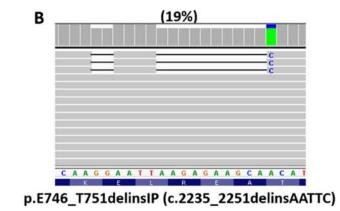

**Figure 1:** *EGFR* mutations detected by next generation sequencing. *EGFR* exon 19 deletion mutation of 3.2% variant allele frequency detected by Integrative Genomics Viewer inspection, but missed by Torrent Variant Caller (**A**). Complex exon 19 deletion mutation composed of two deletions (c.2235 2236del and c.2241 2250del) and one single nucleotide change (c.2251A>C) (**B**).

cDNA change	Amino acid change	Number of cases
Simple deletion		
c.2235_2249del	p.E746_A750del	35
c.2236_2250del	p.E746_A750del	16
c.2240_2254del	p.L747_T751del	3
c.2239_2256del	p.L747_S752del	2
c.2240_2257del	p.L747_P753delinsS	7
Complex deletion		
c.2236_2248delinsGCAC (c.2236_2244del, c.2246A>C, c.2248G>C) <sup>a,b</sup>	p.E746_A750delinsAP	1
c.2235_2251delinsAATTC (c.2235_2236del, c.2241_2250del, c.2251A>C) <sup>b</sup>	p.E746_T751delinsIP	1
c.2237_2252delinsT (c.2237_2251del, c.2252C>T)b	p.E746_T751delinsV	1
c.2237_2255delinsT (c.2237_2254del, c.2255C>T)	p.E746_S752delinsV	2
c.2237_2257delinsTCT (c.2237_2253del, c.2257del) <sup>b</sup>	p.E746_P753delinsVS	1
c.2239_2248delinsC (c.2236_2244del, c.2248G>C)	p.L747_A750delinsP	3
c.2238_2248delinsGC (c.2236_2244del, c.2247A>G, c.2248G>C)	p.L747_A750delinsP	1
c.2239_2251delinsC (c.2238_2249del, c.2251A>C)	p.L747_T751delinsP	2
c.2239_2252delinsCA (c.2239_2248del, c.2251_2252del) <sup>a,b</sup>	p.L747_T751delinsQ	1
c.2239_2258delinsCA (c.2239_2256del, c.2258C>A)	p.L747_P753delinsQ	1
c.2240_2264>CGAAAGG (c.2240_2257del, c.2264C>G) <sup>a,b</sup>	p.L747_A755delinsSKG	1

Table 3: Simple and complex exon 19 deletion mutations of the EGFR gene

<sup>a</sup>Undetectable or detectable with a lower analytic sensitivity by cobas EGFR test.

<sup>b</sup>Undetectable or detectable with a lower analytic sensitivity by therascreen EGFR test.

mutations (Table 2). VAFs were 5% or less in 2.2% of *EGFR* mutations and 10% or less in 7.8% of *EGFR* mutations. VAFs of the *KRAS* gene (27% ± 18%, mean ± standard deviation, P < 0.001 by Student's *t*-test) were significantly lower than those of the *EGFR* gene (35% ± 20%) and significantly higher than those of the *BRAF* (20% ± 16%, P < 0.001) and *PIK3CA* genes (20% ± 19%, P = 0.005). The proportions of *BRAF*-mutated specimens and *PIK3CA*-mutated specimens with < 10% or < 20% VAFs were significantly higher than those of *KRAS*-mutated specimens and *EGFR*-mutated specimens (Table 2).

# Doublet (compound) EGFR mutations

There were 29 tumors with two *EGFR* mutations including 20 (69%) tumors with two non-p.L858R missense mutations, 7 (24%) with p.L858R and non-p. L858R missense mutation, one with p.L747\_P753delinsS and p.T790M and one with p.L747\_T751del and p.K754Q (Table 4). VAFs of doublet mutations were highly concordant (Figure 2A, r = 0.80) except for a tumor with 7.6% p.L858R and 65% p.T790M in a context of 11–30% estimated tumor cellularity, suggesting a germline p.T790M mutation. Doublet mutations were uncommon in tumors with exon 19 deletion (2 of 78, 2.6%) and exon 20 insertion (0 of 18) (Figure 3). Among the single nucleotide mutations, p.L858R (12%) showed a significantly lower incidence of concomitant EGFR mutations than other single nucleotide mutations, 100% in codons 289, 709, 768 and 790 (all *P* values < 0.001) and 70% to 75% in codons 108 (P = 0.01), 719 (P < 0.001) and 861 (P = 0.002) (Figure 3). These included 6 tumors with p.S768I and p.G719A/C/S, one with p.S768I and p.L858R and 4 with the resistance mutation, p.T790M (two with p.L858R, one with p.G719A, and one with exon 19 deletion). There were no coexisting mutations between codon 719 mutations, exon 19 deletion, exon 20 insertion, and p.L858R. All the partners of codon 719 mutations were those less common missense mutations, including 3 tumors with codon 709 mutations of the same alleles. Doublet mutations were observed in 7 of 8 tumors with mutations involving the extracellular domain (codons 108 and 289) including 3 tumors with p.L858R, two with codon 724 mutation, one with p.L861Q, and one with p.G719C.

# Concomitant EGFR and PIK3CA mutations


There were 7 tumors with both *EGFR* and *PIK3CA* mutations. In contrast to double *EGFR* mutations, *EGFR* 

VAFs were not concordant with those of the coexisting *PIK3CA* mutations (Figure 2B, r = 0.07). VAFs of the *EGFR* mutations were equivalent to or higher than those of the coexisting *PIK3CA* mutations (P = 0.009 by paired Student's *t*-test).

#### Concomitant BRAF and KRAS mutations

Activating *KRAS* mutation was detected in 11 (18%) of 63 *BRAF*-mutated tumors, including one (3.0%)

of 33 tumors with a kinase-activated *BRAF* mutation, 5 (31%) of 16 with a kinase-impaired *BRAF* mutation, and 5 (36%) of 14 with a *BRAF* mutation of unknown kinase activity (Table 5) [19–22]. None of the tumors with concomitant *BRAF* and *KRAS* mutations had the p.V600E mutation. The incidence of concomitant *BRAF* mutations with activating *KRAS* mutations is significantly lower in tumors with kinase-activated mutants than those with kinase-impaired (P = 0.01) or kinase-unknown mutants (P = 0.006).

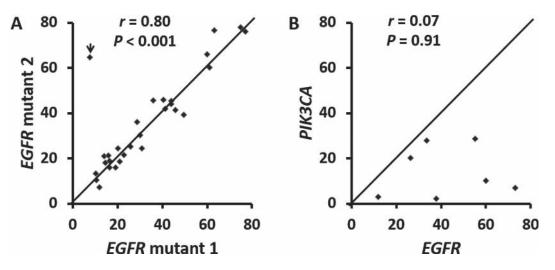

**Figure 2:** Correlation of variant allele frequencies in 29 tumors with doublet (compound) *EGFR* mutations (**A**) and 7 tumors with concomitant *EGFR* and *PIK3CA* mutations (**B**). The r is 0.94 when an outlier (arrowhead in A) was removed. This outlier represents a tumor with a p.T790M germline mutation and a p.L858 somatic mutation. r: Spearman's rank correlation coefficiency.

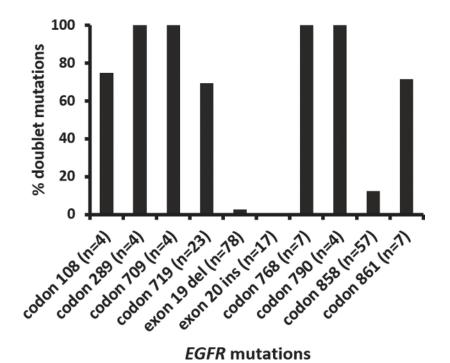


Figure 3: Incidence of doublet (compound) *EGFR* mutations involving different codons. Exon 19 del: exon 19 deletion mutations. Exon 20 insertion/duplication mutations.

mut 1 <sup>a</sup>	Exon (mut 1)	mut 2 <sup>a</sup>	Exon (mut 2)
p.G724A (17%)	18	p.R108G (16%)	3
p.G719C (44%)	18	p.A289V (45%)	7
p.G724S (50%)	18	p.A289V (39%)	7
p.G719A (23%) <sup>b</sup>	18	р.І706Т (22%) <sup>ь</sup>	18
p.G719A (44%) <sup>b</sup>	18	p.E709A (44%) <sup>b</sup>	18
p.G719A (26%) <sup>b</sup>	18	p.E709K (25%) <sup>b</sup>	18
p.G719C (61%) <sup>b</sup>	18	p.E709A (60%) <sup>b</sup>	18
p.G719A (36%)	18	p.S768I (46%)	20
p.G719A (63%)	18	p.G779C (77%)	20
p.G719A (40%)	18	p.S768I (46%)	20
p.G719A (19%)	18	p.T790M (16%)	20
p.G719C (77%)	18	p.S768I (76%)	20
p.G719C (41%)	18	p.S768I (42%)	20
p.G719S (11%)	18	p.S768I (10%)	20
p.G719S (46%)	18	p.S768I (42%)	20
p.G719A (60%)	18	p.L861Q (66%)	21
p.G719D (21%)	18	p.L861Q (19%)	21
p.G719S (31%)	18	p.L861Q (25%)	21
p.L747_T751del (30%) <sup>b,c</sup>	19	p.K754Q (30%) <sup>b,c</sup>	19
p.L747_P753delinsS (12%)	19	p.T790M (7.3%)	20
p.L858R (16%)	21	p.R108K (21%)	3
p.L858R (14%)	21	p.R108K (21%)	3
p.L858R (15%)	21	p.A289V (18%)	7
p.L861Q (75%)	21	p.A289T (78%)	7
p.L858R (10%)	21	p.E709K (13%)	18
p.L858R (29%)	21	p.I744M (36%)	19
p.L858R (7.6%)	21	p.T790M (65%)	20
p.L858R (20%)	21	p.S768I (24%)	20
p.L861Q (17%)	21	p.T790M (19%)	20

Table 4: Doublet (compound) EGFR mutations in 1006 lung cancers

<sup>a</sup>Percentage in the parenthesis indicates variant allele frequency.

<sup>b</sup>Located within the same allele.

<sup>c</sup>Double *EGFR* mutations seen in this tumor may also be categorized as a complex exon 19 deletion (c.2238\_2260delinsATCTCCGC, p.L747\_K754delinsSPQ).

# DISCUSSION

NGS has been clinically validated for mutational profiling of lung cancer [7, 23, 24]. We have previously shown a test feasibility of 94% among the first 625 lung cancer specimens submitted for NGS testing. This included an approximately 3% of specimens rejected due to inadequate specimens and approximately 3% of specimens failed the NGS assay [14]. This retrospective analysis of 1006 lung cancers reaffirms the strength of NGS in clinical mutational profiling. NGS demonstrates a great analytic sensitivity, broad reportable ranges,

and accurate detection and annotation of complex indel mutations. With an analytic sensitivity of 10-20% VAF, Sanger sequencing would have missed 7.8% or 27% of *EGFR* mutations with less than 10% or 20% VAFs in this series. The analytic sensitivity can be improved to approximately 1-5% VAF by mutation-specific real time PCR assays such as cobas EGFR mutation test and therascreen EGFR RGQ PCR Kit, which were designed to detect only hot spot *EGFR* mutations [9–13].

NGS detected a variety of uncommon mutations located outside the reportable ranges of cobas and therascreen tests, including 4 codon 108 mutations

	8	0
mut 1 (VAF%)	mut 2 (VAF%)	<b>BRAF</b> kinase activity
<i>KRAS</i> p.Q61H (12%)	BRAF p.G464V (9.5%)	activated <sup>a</sup> [19]
<i>KRAS</i> p.G12D (19%)	BRAF p.G466E (3.3%)	Impaired [19]
<i>KRAS</i> p.A146T (30%)	<i>BRAF</i> p.G466V (35%)	Impaired [19]
<i>KRAS</i> p.G12D (19%)	<i>BRAF</i> p.G466V (14%)	Impaired [19]
<i>KRAS</i> p.G12V (14%)	BRAF p.D594G (2.6%)	Impaired [20]
<i>KRAS</i> p.Q61H (24%) <sup>b</sup>	<i>BRAF</i> p.D594N (5.7%) <sup>b</sup>	Impaired [21]
<i>KRAS</i> p.G12V (16%) <sup>b</sup>	<i>BRAF</i> p.T440I (13%) <sup>b</sup>	unknown
KRAS p.G12C (14%)	BRAF p.G466L (12%)	unknown
<i>KRAS</i> p.G12R (12%) <sup>b</sup>	<i>BRAF</i> p.G469R (11%) <sup>b</sup>	unknown
<i>KRAS</i> p.G13D (25%) <sup>b</sup>	<i>BRAF</i> p.G469V (20%) <sup>b</sup>	unknown
<i>KRAS</i> p. A146T (11%)	BRAF p.K601N (11%)	unknown

Table 5: Concomitant BRAF mutation and activating KRAS mutation in lung cancers

VAF: variant allele frequency.

<sup>a</sup>Categorized as intermediate activity mutants by Wan et al. [19].

<sup>b</sup>These 4 cases have been reported previously [22].

and 4 codon 289 mutations involving the extracellular domain of EGFR. These mutations are more prevalent in glioblastomas and their clinical significance in lung cancer is not known [25]. NGS also detected 3 complex exon 19 mutations not included in the design of cobas and therascreen tests and additional 3 not included in the therascreen test. Among 18 tumors with TKI-resistant exon 20 insertion mutations detected by NGS, cobas test would have detected only 8 with c.2300\_2308dup or c.2317\_2319dup, and therascreen test would have detected only 3 with c.2317\_2319dup. Both cobas and therascreen tests would also have missed TKI-sensitive exon 20 insertion mutations, such as p.A763\_Y764insFQEA, which can be detected by this NGS assay as shown previously [26, 27].

Complex exon 19 deletion mutations may not be accurately characterized by Sanger sequencing [28], partly because the accompanied single nucleotide change or second deletion change may be difficult to interpret in the presence of an underlying deletion mutation, especially at a lower level of VAF. Furthermore, Sanger sequencing cannot distinguish if the two adjacent sequence changes are located within the same allele or different alleles without laborious cloning of the PCR amplicons for sequencing. NGS platforms provide individual sequencing information from a single molecule and, therefore, can confirm that the two sequence changes are always located within the same allele to form a complex exon 19 deletion. With accurate detection and annotation of the complex exon 19 deletion, further studies are warranted to elucidate if the point mutation component of the complex exon 19 deletion mutations may decrease TKI efficacy.

Application of assays with broader reportable ranges may shed light on the significance of uncommon

mutations. For example, mutations involving codon 33 of the KRAS gene were recently found to be oncogenic [29]. NGS detected 17 unique *PIK3CA* mutations including 3 novel ones. Commonly reported mutations involving codons 542, 545 and 1047 accounted for only 65% of the PIK3CA mutations. A total of 24 unique BRAF mutations including 4 novel ones were detected among 63 BRAFmutated lung cancers. The p.V600E mutant accounted for only 24% of BRAF mutations while kinase-impaired BRAF mutants involving codons 466 and 594 were seen in 25% of BRAF mutations. In a previous study of combined lung cancer, colorectal cancer and melanoma specimens, kinase-impaired BRAF mutants were associated with a higher incidence of a concomitant activating KRAS/ NRAS mutation [22]. This is confirmed by a larger cohort of lung cancer specimens in this study. In vitro studies have shown that in the presence of oncogenic RAS proteins, kinase-impaired BRAF forms a complex with CRAF and leads to hyperactivation of the CRAF/MEK/ ERK cascade, suggesting MEK inhibitors or CRAF inhibitors may benefit patients with concomitant kinaseimpaired BRAF mutation and activating RAS mutation [20, 21]. Dabrafenib (a selective BRAF inhibitor) alone or combined with trametinib (a MEK inhibitor) has shown efficacy in p.V600E-mutated lung cancers [30, 31]. The NCI-Molecular Analysis for Therapy Choice (NCI-MATCH) Trial also includes an arm for targeting tumors with non-p.V600E BRAF mutations with trametinib (https://www.cancer.gov/about-cancer/treatment/clinicaltrials/nci-supported/nci-match, accessed 1/19/2017).

False negative results were a major concern when NGS platforms were initially implemented in the clinical diagnostic setting. During our clinical validation of this NGS platform, we have found that Torrent Variant Caller alone may miss EGFR p.L858R [18]. Therefore, we have included IGV inspection of each amplicon as a second analysis pipeline. We also examined the indel mutations within EGFR exons 19 and 20 by a sizing assay. All EGFR mutations, missed by older versions of Torrent Variant Caller, were observed by IGV. Although recent versions did not miss mutations in lung cancers, version 5.0.2.1 (since December 2015) did miss a 8.4% KIT p.K558 D572del (45-base deletion mutation) in a gastrointestinal stromal tumor specimen and version 5.0.4.0 (since September 2016) missed a 2.6% BRAF p.V600E in a melanoma specimen with less than 10% estimated tumor cellularity (unpublished data). All exons 19 and 20 indel mutations detected by the sizing assay were also observed by IGV inspection. The results indicate that indel mutations of 18 or less bases can be reliably detected by Torrent Suite analysis combined with direct visual inspection of the binary sequence alignment/map file using IGV. Longer indel mutations such as internal tandem duplication mutations of the FLT3 gene may not be detected by NGS assays without bioinformatics pipelines designed for longer indel mutations [32, 33].

Consistent VAFs over a 4-year period in positive control specimens highlighted the precise quantitative nature of NGS assays [27, 34]. Analysis of VAF may vield important information regarding mutant allelespecific imbalance (such as gene amplification or loss of heterozygosity), tumor heterogeneity, and germline mutations [7, 34, 35]. We have shown that lower than expected VAF indicated tumor heterogeneity while higher than expected VAF indicated mutant allele-specific imbalance [34, 35]. In this study, we found an equivalent higher EGFR VAF than the coexisting PIK3CA VAF. The results suggest a higher incidence of mutant allele-specific imbalance (most likely duplication or amplification) of the EGFR gene or the presence of PIK3CA mutation in a subpopulation of the EGFR-mutated tumors which may contribute to TKI resistance [27]. This was confirmed by subarea analysis of a specimen containing 2.8% PIK3CA mutation and 12% EGFR p.L858R mutation. Only one of 8 fragments showed both mutations with concordant VAFs (data not shown). In contrast, VAFs of the doublet EGFR mutations were highly concordant except for one specimen with 7.6% p.L858R and 65% p.T790M in a context of 11-30% estimated tumor cellularity, suggesting a germ-line p.T790M mutation [36].

Doublet *EGFR* mutations are not uncommon (1.6% to 14%) [37–40]. A summary from 66 publications showed 96 (6%) doublets in 1621 *EGFR*-mutated lung cancers [41], including several tumors which were indeed complex exon 19 deletion. The incidence of doublet *EGFR* mutations is higher by using sensitive assays with broader reportable ranges such as NGS assays. We found two *EGFR* mutations in 29 (16%) of 187 *EGFR*-mutated tumors while 11 tumors with complex exon 19 deletion were excluded. Concordant VAFs between two

mutations suggests one dominant tumor population with two mutations rather than two tumor populations each containing a mutation. Both mutations may be needed to initiate the founder clone or an EGFR-mutated subclone that has become the dominant tumor population after acquiring the second mutation. Concomitant mutations between exon 19 deletion, exon 20 insertion and p.L858R are uncommon [41]. Only one doublet consisting of exon 19 deletion and p.L858R was reported among 1621 EGFR-mutated lung cancers [41]. In contrast, non-p. L858R missense mutations were often seen in doublet mutations, suggesting a lower oncogenic potential of these mutations. In vitro studies have shown a comparable or higher level of catalytic phosphorylating activity in mutants involving conservative codons at 709, 719, 768, 790 and 861 compared to the wild-type, but a lower level of kinase activity with respect to exon 19 deletion or p.L858R mutants [42-45]. Significant enhancement of kinase activity observed in doublet with p.T790M and p.L858R or exon 19 deletion suggests additive oncogenic effect from p.T790M [43, 44].

NGS demonstrates a high analytic sensitivity, quantitative measurement of VAF, single molecule sequencing of complex exon 19 deletion, and broad reportable ranges with simultaneous detection of doublet *EGFR* mutations and concomitant *BRAF* and *KRAS* mutations in the clinical diagnostic setting. Further studies are warranted to elucidate the clinical and/or biological significance of uncommon mutations, doublet non-p. L858R missense mutations of *EGFR*, and concomitant kinase-impaired *BRAF* mutations in lung cancers.

# **MATERIALS AND METHODS**

# Materials

NGS was conducted in 1103 formalin-fixed paraffin-embedded (FFPE) specimens with a diagnosis of adenocarcinoma in situ (5 specimens), invasive adenocarcinoma (1033 specimens), adenosquamous carcinoma (10 specimens) or non-small cell carcinoma (55 specimens) of lung submitted to the Molecular Diagnostics Laboratory at The Johns Hopkins Hospital between April 2013 and June 2016. Specimens with prior TKI therapy were not included. There were 499 (45%) resection specimens, 341 (31%) biopsy specimens, and 204 (19%) fine needle aspiration specimens, 55 (5.0%) pleural or pericardial effusion specimen, 2 bronchial brushing specimen, one bronchoalveolar lavage specimen and one curetting of bone (Supplementary Table 6). Twenty-nine (2.6%) specimens failed (Supplementary Table 6). The remaining 1074 specimens with successful NGS were submitted from 1006 tumors of 987 patients (Supplementary Table 7). Tissue blocks with adequate tumor cellularity were selected by the pathologists who made the diagnosis. One hematoxylin & eosin (H&E) slide

followed by 5-10 unstained slides and one additional H&E slide were prepared with PCR precautions. The H&E slide was examined and marked by PI (pulmonary pathologist), MTL and/or GZ (molecular pathologists) for subsequent macro-dissection of FFPE neoplastic tissues from 3-10 unstained slides of 5 or 10-micron thick sections. DNA was isolated from the designated area(s) using the Pinpoint DNA Isolation System (Zymo Research, Irvine, CA), followed by further purification via the QIAamp Mini Kit (Qiagen, Valencia, CA) [46]. DNA concentration was measured by Qubit 2.0 Fluorometer (Life Technologies, Carlsbad, California). Specimens with concentration less than 10 ng/µl or bony specimens were concentrated with Amicon Ultra 0.5 ml centrifugal filters with ultracel-30K membrane (Millipore Sigma, Darmstadt, Germany) after April 2015. The Johns Hopkins Medicine institutional review board granted approval to this study.

#### Next generation sequencing (NGS)

NGS was conducted using AmpliSeq Cancer Hotspot Panel (v2) for targeted multi-gene amplification, as described previously [18, 34]. Briefly, we used the Ion AmpliSeq Library Kit 2.0 for library preparation, Ion OneTouch 200 Template Kit v2 DL (replaced by Ion Personal Genome Machine Hi-Q OT2 Kit lately) and Ion OneTouch 2 Instrument for emulsion PCR and template preparation, and the Ion Personal Genome Machine 200 Sequencing Kit (replaced by Ion Personal Genome Machine Hi-Q Sequencing Kit lately) with the Ion 318 Chip and Personal Genome Machine as the sequencing platform (Life Technologies). The DNA input ranged from 1 ng to 30 ng, as measured by Qubit 20 Fluorometer (Life Technologies). Up to 8 specimens were barcoded using Ion Xpress Barcode Adapters (Life Technologies) for each Ion 318 chip. One to three controls (a non-template control, a normal peripheral blood control from a male, and/or positive control specimens) were included in each chip. Positive controls were mixed DNA specimens from several cell lines with known mutations as reported previously [27, 34].

Redundant bioinformatics pipelines are essential for NGS analysis, as a single analysis pipeline may give false negative and false positive results. Direct visual inspection of the binary sequence alignment/map file using the Broad Institute's Integrative Genomics Viewer (IGV) (http:// www.broadinstitute.org/igv/) was implemented during the validation processes of this assay after we found that Torrent Variant Caller missed the most common *EGFR* point mutation (p.L858R) in a lung cancer specimen [18, 47]. In this study, sequencing data were analyzed using Torrent Suite (Life Technologies) as described previously [18]. Mutations were identified and annotated through both Torrent Variant Caller (Life Technologies) and direct visual inspection of the binary sequence alignment/map file using

IGV. All specimens were analyzed for AKT1 (NM 005163), BRAF (NM 004333), EGFR (NM 005228), ERBB2 (NM 004448), KRAS (NM 033360), NRAS (NM 002524) and PIK3CA (NM 006218) genes. During our validation of this NGS assay, a cutoff of background noise at 2% was chosen for single nucleotide variant according to a study of 16 non-neoplastic FFPE tissues [18]. We also developed a statistical model to determine the read depth needed for a given percent tumor cellularity and number of functional genomes. With sufficient DNA input, the limit of detection is dictated by the depth of coverage (or number of sequencing reads). Approximately 150 and 500 reads is needed to detect a heterozygous mutation at a 99% confidence in a specimen with 20% and 10% tumor cellularity, respectively. The reportable ranges and reference ranges for the 7 genes have been reported previously [27, 34]. BRAF mutation data before September 2014 has also been published previously together with colorectal cancer specimens and melanoma specimens [22].

# Alternative assays

Insertion/deletion mutations within exons 19 and 20 of the *EGFR* genes were also examined by sizing using capillary electrophoresis. PCR followed by capillary electrophoresis was conducted as described previously [18]. This was discontinued after April 2014. Mutations not reported in the COSMIC database were confirmed by either Sanger sequencing or pyrosequencing. This policy was also discontinued after April 2015.

# Statistical analysis

Student's *t*-test,  $\chi^2$  test or Fisher exact test was performed to calculate *P* values. Correlation of frequencies between two mutations was examined by Spearman's rank correlation coefficient (denoted as *r*) using the GraphPad Prism software (GraphPad Software, ver5, La Jolla, CA) as described previously [48].

# ACKNOWLEGMENTS

We would like to thank Kevin Zhen for assistance in preparation of figures.

# **CONFLICTS OF INTEREST**

The authors declare no conflicts of interest.

# FUNDING

This retrospective analysis for quality improvement was supported by 1UM1CA186691-01 from the NIH-National Cancer Institute.

# **REFERENCE**S

- Lynch TJ, Bell DW, Sordella R, Gurubhagavatula S, Okimoto RA, Brannigan BW, Harris PL, Haserlat SM, Supko JG, Haluska FG, Louis DN, Christiani DC, Settleman J, Haber DA. Activating mutations in the epidermal growth factor receptor underlying responsiveness of non-small-cell lung cancer to gefitinib. N Engl J Med. 2004; 350:2129–2139.
- Paez JG, Jänne PA, Lee JC, Tracy S, Greulich H, Gabriel S, Herman P, Kaye FJ, Lindeman N, Boggon TJ, Naoki K, Sasaki H, Fujii Y, et al. EGFR mutations in lung cancer: correlation with clinical response to gefitinib therapy. Science. 2004; 304:1497–1500.
- Sequist LV, Yang JC, Yamamoto N, O'Byrne K, Hirsh V, Mok T, Geater SL, Orlov S, Tsai CM, Boyer M, Su WC, Bennouna J, Kato T, et al. Phase III study of afatinib or cisplatin plus pemetrexed in patients with metastatic lung adenocarcinoma with EGFR mutations. J Clin Oncol. 2013; 31:3327–3334.
- 4. Keedy VL, Temin S, Somerfield MR, Beasley MB, Johnson DH, McShane LM, Milton DT, Strawn JR, Wakelee HA, Giaccone G. American Society of Clinical Oncology provisional clinical opinion: epidermal growth factor receptor (EGFR) Mutation testing for patients with advanced non-small-cell lung cancer considering first-line EGFR tyrosine kinase inhibitor therapy. J Clin Oncol. 2011; 29:2121–2127.
- 5. Lindeman NI, Cagle PT, Beasley MB, Chitale DA, Dacic S, Giaccone G, Jenkins RB, Kwiatkowski DJ, Saldivar JS, Squire J, Thunnissen E, Ladanyi M. Molecular testing guideline for selection of lung cancer patients for EGFR and ALK tyrosine kinase inhibitors: guideline from the College of American Pathologists, International Association for the Study of Lung Cancer, and Association for Molecular Pathology. Arch Pathol Lab Med. 2013; 137: 828–860.
- Ellison G, Zhu G, Moulis A, Dearden S, Speake G, McCormack R. EGFR mutation testing in lung cancer: a review of available methods and their use for analysis of tumour tissue and cytology samples. J Clin Pathol. 2013; 66:79–89.
- Frampton GM, Fichtenholtz A, Otto GA, Wang K, Downing SR, He J, Schnall-Levin M, White J, Sanford EM, An P, Sun J, Juhn F, Brennan K, et al. Development and validation of a clinical cancer genomic profiling test based on massively parallel DNA sequencing. Nat Biotechnol. 2013; 31:1023–1031.
- Chen G, Yang Z, Eshleman JR, Netto GJ, Lin MT. Molecular diagnostics for precision medicine in colorectal cancer: current status and future perspective. Biomed Res Int. 2016;2016:9850690.
- O'Donnell P, Ferguson J, Shyu J, Current R, Rehage T, Tsai J, Christensen M, Tran HB, Chien SS, Shieh F, Wei W, Lawrence HJ, Wu L, et al. Analytic performance studies and clinical reproducibility of a real-time PCR assay for the

detection of epidermal growth factor receptor gene mutations in formalin-fixed paraffin-embedded tissue specimens of non-small cell lung cancer. BMC Cancer. 2013; 13:210.

- Kimura H, Ohira T, Uchida O, Matsubayashi J, Shimizu S, Nagao T, Ikeda N, Nishio K. Analytical performance of the cobas EGFR mutation assay for Japanese non-small-cell lung cancer. Lung Cancer. 2014; 83:329–333.
- Wong AT, To RM, Wong CL, Chan WK, Ma ES. Evaluation of 2 real-time PCR assays for *in vitro* diagnostic use in the rapid and multiplex detection of EGFR gene mutations in NSCLC. Diagn Mol Pathol. 2013; 22:138–143.
- 12. Vallée A, Le Loupp AG, Denis MG. Efficiency of the Therascreen<sup>®</sup> RGQ PCR kit for the detection of EGFR mutations in non-small cell lung carcinomas. Clin Chim Acta. 2014; 429:8–11.
- Syed YY. therascreen<sup>®</sup> EGFR RGQ PCR Kit: A companion diagnostic for afatinib and gefitinib in non-small cell lung cancer. Mol Diagn Ther. 2016; 20:191–198.
- Zheng G, Tsai H, Tseng LH, Illei P, Gocke CD, Eshleman JR, Netto G, Lin MT. Test feasibility of nextgeneration sequencing assays in clinical mutation detection of small biopsy and fine needle aspiration specimens. Am J Clin Pathol. 2016; 145:696–702.
- Sherwood JL, Müller S, Orr MC, Ratcliffe MJ, Walker J. Panel based MALDI-TOF tumour profiling is a sensitive method for detecting mutations in clinical non small cell lung cancer tumour. PLoS One. 2014; 9:e100566.
- 16. Quinn AM, Hickson N, Adaway M, Priest L, Jaeger E, Udar N, Keeling C, Kamieniorz M, Dive C, Wallace A, Byers RJ, Newman WG, Nonaka D, Blackhall FH. Diagnostic mutation profiling and validation of nonsmall-cell lung cancer small biopsy samples using a high throughput platform. J Thorac Oncol. 2015; 10:784–792.
- 17. Dienstmann R, Rodon J, Barretina J, Tabernero J. Genomic medicine frontier in human solid tumors: prospects and challenges. J Clin Oncol. 2013; 31:1874–1884.
- Lin MT, Mosier SL, Thiess M, Beierl KF, Debeljak M, Tseng LH, Chen G, Yegnasubramanian S, Ho H, Cope L, Wheelan SJ, Gocke CD, Eshleman JR. Clinical validation of KRAS, BRAF, and EGFR mutation detection using nextgeneration sequencing. Am J Clin Pathol. 2014; 141:856–866.
- Wan PT, Garnett MJ, Roe SM, Lee S, Niculescu-Duvaz D, Good VM, Jones CM, Marshall CJ, Springer CJ, Barford D, Marais R; Cancer Genome Project. Mechanism of activation of the RAF-ERK signaling pathway by oncogenic mutations of B-RAF. Cell. 2004; 116:855–867.
- Smalley KS, Xiao M, Villanueva J, Nguyen TK, Flaherty KT, Letrero R, Van Belle P, Elder DE, Wang Y, Nathanson KL, Herlyn M. CRAF inhibition induces apoptosis in melanoma cells with non-V600E BRAF mutations. Oncogene. 2009; 28:85–94.
- Heidorn SJ, Milagre C, Whittaker S, Nourry A, Niculescu-Duvas I, Dhomen N, Hussain J, Reis-Filho JS, Springer CJ, Pritchard C, Marais R. Kinase-dead BRAF and oncogenic

RAS cooperate to drive tumor progression through CRAF. Cell. 2010; 140:209–221.

- Zheng G, Tseng LH, Chen G, Haley L, Illei P, Gocke CD, Eshleman JR, Lin MT. Clinical detection and categorization of uncommon and concomitant mutations involving BRAF. BMC Cancer. 2015;15:779.
- de Biase D, Visani M, Malapelle U, Simonato F, Cesari V, Bellevicine C, Pession A, Troncone G, Fassina A, Tallini G. Next-generation sequencing of lung cancer EGFR exons 18–21 allows effective molecular diagnosis of small routine samples (cytology and biopsy). PLoS One. 2013; 8:e83607.
- 24. Tafe LJ, Pierce KJ, Peterson JD, de Abreu F, Memoli VA, Black CC, Pettus JR, Marotti JD, Gutmann EJ, Liu X, Shirai K, Dragnev KH, Amos CI, Tsongalis GJ. Clinical genotyping of non-small cell lung cancers using targeted next-generation sequencing: utility of identifying rare and co-mutations in oncogenic driver genes. Neoplasia. 2016; 18:577–583
- 25. Lee JC, Vivanco I, Beroukhim R, Huang JH, Feng WL, DeBiasi RM, Yoshimoto K, King JC, Nghiemphu P, Yuza Y, Xu Q, Greulich H, Thomas RK, et al. Epidermal growth factor receptor activation in glioblastoma through novel missense mutations in the extracellular domain. PLoS Med. 2006; 3:e485.
- 26. Yasuda H, Park E, Yun CH, Sng NJ, Lucena-Araujo AR, Yeo WL, Huberman MS, Cohen DW, Nakayama S, Ishioka K, Yamaguchi N, Hanna M, Oxnard GR, et al. Structural, biochemical, and clinical characterization of epidermal growth factor receptor (EGFR) exon 20 insertion mutations in lung cancer. Sci Transl Med. 2013; 5:216ra177.
- Belchis DA, Tseng LH, Gniadek T, Haley L, Lokhandwala P, Illei P, Gocke CD, Forde P, Brahmer J, Askin FB, Eshleman JR, Lin MT. Heterogeneity of resistance mutations detectable by next-generation sequencing in TKI-treated lung adenocarcinoma. Oncotarget. 2016; 7:45237–45248. https://doi.org/10.18632/oncotarget.9931.
- 28. Marchetti A, Del Grammastro M, Filice G, Felicioni L, Rossi G, Graziano P, Sartori G, Leone A, Malatesta S, Iacono M, Guetti L, Viola P, Mucilli F, et al. Complex mutations & subpopulations of deletions at exon 19 of EGFR in NSCLC revealed by next generation sequencing: potential clinical implications. PLoS One. 2012; 7:e42164.
- 29. Kim E, Ilic N, Shrestha Y, Zou L, Kamburov A, Zhu C, Yang X, Lubonja R, Tran N, Nguyen C, Lawrence MS, Piccioni F, Bagul M, et al. Systematic functional interrogation of rare cancer variants identifies oncogenic alleles. Cancer Discov. 2016; 6:714–726.
- 30. Planchard D, Kim TM, Mazieres J, Quoix E, Riely G, Barlesi F, Souquet PJ, Smit EF, Groen HJ, Kelly RJ, Cho BC, Socinski MA, Pandite L, et al. Dabrafenib in patients with BRAF(V600E)-positive advanced non-small-cell lung cancer: a single-arm, multicentre, open-label, phase 2 trial. Lancet Oncol. 2016; 17:642–650.
- 31. Planchard D, Besse B, Groen HJ, Souquet PJ, Quoix E, Baik CS, Barlesi F, Kim TM, Mazieres J, Novello S,

Rigas JR, Upalawanna A, D'Amelio AM Jr, et al. Dabrafenib plus trametinib in patients with previously treated BRAF(V600E)-mutant metastatic non-small cell lung cancer: an open-label, multicentre phase 2 trial. Lancet Oncol. 2016; 17:984–993.

- Spencer DH, Abel HJ, Lockwood CM, Payton JE, Szankasi P, Kelley TW, Kulkarni S, Pfeifer JD, Duncavage EJ. Detection of FLT3 internal tandem duplication in targeted, short-read-length, next-generation sequencing data. J Mol Diagn. 2013; 15:81–93.
- 33. Kadri S, Zhen CJ, Wurst MN, Long BC, Jiang ZF, Wang YL, Furtado LV, Segal JP. Amplicon indel hunter is a novel bioinformatics tool to detect large somatic insertion/ deletion mutations in amplicon-based next-generation sequencing data. J Mol Diagn. 2015; 17:635–643.
- Haley L, Tseng LH, Zheng G, Dudley J, Anderson DA, Azad NS, Gocke CD, Eshleman JR, Lin MT. Performance characteristics of next-generation sequencing in clinical mutation detection of colorectal cancers. Mod Pathol. 2015; 28:1390–1399.
- 35. Dudley J, Tseng LH, Rooper L, Harris M, Haley L, Chen G, Gocke CD, Eshleman JR, Lin MT. Challenges posed to pathologists in the detection of KRAS mutations in colorectal cancers. Arch Pathol Lab Med. 2015; 139:211–218.
- 36. Bell DW, Gore I, Okimoto RA, Godin-Heymann N, Sordella R, Mulloy R, Sharma SV, Brannigan BW, Mohapatra G, Settleman J, Haber DA. Inherited susceptibility to lung cancer may be associated with the T790M drug resistance mutation in EGFR. Nat Genet. 2005; 37:1315–1316.
- Wu SG, Chang YL, Hsu YC, Wu JY, Yang CH, Yu CJ, Tsai MF, Shih JY, Yang PC. Good response to gefitinib in lung adenocarcinoma of complex epidermal growth factor receptor (EGFR) mutations with the classical mutation pattern. Oncologist. 2008; 13:1276–1284.
- Kobayashi S, Canepa HM, Bailey AS, Nakayama S, Yamaguchi N, Goldstein MA, Huberman MS, Costa DB. Compound EGFR mutations and response to EGFR tyrosine kinase inhibitors. J Thorac Oncol. 2013; 8:45–51.
- 39. Li S, Li L, Zhu Y, Huang C, Qin Y, Liu H, Ren-Heidenreich L, Shi B, Ren H, Chu X, Kang J, Wang W, Xu J, et al. Coexistence of EGFR with KRAS, or BRAF, or PIK3CA somatic mutations in lung cancer: a comprehensive mutation profiling from 5125 Chinese cohorts. Br J Cancer. 2014; 110:2812–2820.
- Peng L, Song Z, Jiao S. Comparison of uncommon EGFR exon 21 L858R compound mutations with single mutation. Onco Targets Ther. 2015; 8:905–910.
- 41. Chen Z, Feng J, Saldivar JS, Gu D, Bockholt A, Sommer SS. EGFR somatic doublets in lung cancer are frequent and generally arise from a pair of driver mutations uncommonly seen as singlet mutations: one-third of doublets occur at five pairs of amino acids. Oncogene. 2008; 27:4336–4343.
- 42. Chen YR, Fu YN, Lin CH, Yang ST, Hu SF, Chen YT, Tsai SF, Huang SF. Distinctive activation patterns in

constitutively active and gefitinib-sensitive EGFR mutants. Oncogene. 2006; 25:1205–1215.

- 43. Mulloy R, Ferrand A, Kim Y, Sordella R, Bell DW, Haber DA, Anderson KS, Settleman J. Epidermal growth factor receptor mutants from human lung cancers exhibit enhanced catalytic activity and increased sensitivity to gefitinib. Cancer Res. 2007; 67:2325–2330.
- 44. Godin-Heymann N, Bryant I, Rivera MN, Ulkus L, Bell DW, Riese DJ 2nd, Settleman J, Haber DA. Oncogenic activity of epidermal growth factor receptor kinase mutant alleles is enhanced by the T790M drug resistance mutation. Cancer Res. 2007; 67:7319–7326.
- 45. Yun CH, Boggon TJ, Li Y, Woo MS, Greulich H, Meyerson M, Eck MJ. Structures of lung cancer-derived EGFR mutants and inhibitor complexes: mechanism of activation and insights into differential inhibitor sensitivity. Cancer Cell. 2007; 11:217–227.

- 46. Lin MT, Tseng LH, Rich RG, Hafez MJ, Harada S, Murphy KM, Eshleman JR, Gocke CD. Δ-PCR, A simple method to detect translocations and insertion/deletion mutations. J Mol Diagn. 2011; 13:85–92.
- Thorvaldsdóttir H, Robinson JT, Mesirov JP. Integrative Genomics Viewer (IGV): high-performance genomics data visualization and exploration. Brief Bioinform. 2013; 14:178–192.
- 48. Dudley JC, Gurda GT, Tseng LH, Anderson DA, Chen G, Taube JM, Gocke CD, Eshleman JR, Lin MT. Tumor cellularity as a quality assurance measure for accurate clinical detection of BRAF mutations in melanoma. Mol Diagn Ther. 2014; 18:409–418.