

Comparison of US Pharmacopeia Simulated Intestinal Fluid TS (without pancreatin) and Phosphate Standard Buffer pH 6.8, TS of the International Pharmacopoeia with Respect to Their Use in *In Vitro* Dissolution Testing

email:dressman@em.uni-frankfurt.de

Erika Stippler^{1,2}, Sabine Kopp^{3,4}, Jennifer B. Dressman^{1,5}

Introduction

Buffer solutions suitable for dissolution tests are described as such only in The United States Pharmacopeia (USP); The International Pharmacopoeia (IntPh), for example, broadly describes the composition of buffer solutions without specifying their fields of applicability. In this study Simulated Intestinal Fluid (SIFsp), which is described in The 26th United States Pharmacopeia (USP 26) as an 0.05 M buffer solution containing potassium dihydrogen phosphate (the pancreatin is omitted from the composition for most dissolution testing applications), and Phosphate Standard Buffer pH 6.8 TS, which is described in Volume 1 of The International Pharmacopoeia, Third Edition, (IntPh 3) as a mixture of 0.025 M potassium dihydrogen phosphate buffer solution and 0.025 M di-sodium hydrogen phosphate buffer solution, were investigated with respect to their interchangeability as dissolution test media.

Both are pH 6.8 phosphate buffers, and they have similar osmolalities (114 mOsmol/kg) and buffer capacities (18 mEq/L/pH/unit). However, in the USP medium potassium (50mM) dominates over sodium (22mM) whereas in the IntPh buffer sodium (49.7mM) dominates over potassium (25mM). [Note: For the current investigations the IntPh phosphate buffer media was used as monographed in the 3rd edition, Volume 1. With the publication of Volume 5 the composition of the buffer pH 6.8 medium has been changed. Volume 5 describes a 0.04 M disodium hydrogen phosphate buffer solution for Buffer pH 6.8.].

The applicability and the interchangeability of these two compositionally similar but not identical buffers as dissolution media was investigated for products containing one of three drug substances: Metronidazole, a typical BCS Class 1 drug, and indometacin and ibuprofen, two drugs belonging to

BCS Class 2. In each case several commercial products currently available on the German market were investigated. Products included uncoated, film coated, and sugar coated tablets as well as hard gelatin capsules.

Experimental Investigations

Materials

Chemical reference standard substances (CRS)

Metronidazole CRS (lot 00616) was purchased from Aldrich (Germany), indometacin CRS (lot 0111A079) was purchased from Synopharm (Germany), and ibuprofen CRS (lot 10337012) was purchased from Caelo (Germany).

Reagents

All reagents used were of analytical grade and were purchased from VWR (Germany).

Dissolution media

The compositions of dissolution media used for in vitro dissolution testing are described in Table 1. The reagents were added to the deionized water under stirring to yield clear solutions. The pH value of the dissolution media was checked and adjusted as necessary to $\text{pH } 6.8 \pm 0.05$ with 0.1 M sodium hydroxide or 0.1 M phosphoric acid solution.

Pharmaceutical preparations

The pharmaceutical preparations (all conventional immediate release solid oral dosage forms) investigated are described in Table 2.

Analytical instruments

Dissolution tests were performed with the USP 26 apparatus 2, using an Erweka model DT6R (Erweka, Germany) apparatus. For manual sampling 5 mL syringes fitted with stainless tubing were used. Samples were filtered immediately after sampling through PTFE filter sticks with 45 μm

¹J. W. Goethe University Frankfurt/Main, Department of Pharmaceutical Technology, Germany

²PHAST GmbH, Homburg/Saar, Germany

³World Health Organization, Geneva, Switzerland

⁴The opinions expressed in this article are those of the author and

do not necessarily represent the decisions or stated policy of the World Health Organization.

⁵Corresponding author, J. W. Goethe University Frankfurt/Main, Department of Pharmaceutical Technology, Biozentrum, Marie-Curie-Str. 9, 60439 Frankfurt, Germany

Table 1: Composition of dissolution media for *in vitro* dissolution testing

Medium	Composition	
Simulated Intestinal Fluid, pH 6.8 ¹ (SIFsp); USP 26	KH ₂ PO ₄	68.05 g
	NaOH	8.96 g
	deionized water to	10.0 L
Phosphate Standard Buffer pH 6.8 ² IntPh 3	KH ₂ PO ₄	34.0 g
	Na ₂ HPO ₄	35.3 g
	deionized water to	10.0 L

¹ The composition of SIFsp is the same as in USP 25.
² The composition of the Phosphate Standard Buffer pH 6.8 has been changed with the publication of Volume 5 of the IntPh.

Table 2: Pharmaceutical preparations under investigation

Drug substance	Drug product	Technol. description	Expiry date
Ibuprofen	IBU-Product A	film coated tablets	11.2002
	IBU-Product B	sugar coated tablets	08.2004
	IBU-Product C	capsules	09.2004
Indometacin	INDO-Product A	uncoated tablets	07.2003
	INDO-Product B	capsules	11.2004
Metronidazole	METRO-Product A	film coated tablets	06.2005
	METRO-Product B	uncoated tablets	09.2004

Table 3: Description of the methods for dissolution testing

Parameter of the method	Method 1	Method 2
Apparatus	paddle	paddle
Dissolution medium	Simulated Intestinal Fluid, pH 6.8 (USP 26)	Phosphate Standard Buffer pH 6.8 (IntPh 3)
Degassing	yes	yes
Volume	500 mL	500 mL
Agitation	75 rpm	75 rpm
Temperature	37°C	37°C
Sampling time	10, 20, 30, 45, 60 min	10, 20, 30, 45, 60 min

pore size. The dissolution system used was calibrated to meet USP apparatus suitability test requirements.

A UV-Vis spectrophotometer model DU 640i (Beckmann Coulter, Germany) was used to analyze the samples for dissolved drug. The instrument was calibrated periodically according to the Pharmacopeial requirements.

For the pH measurements of dissolution media an MP 230 pH meter (Mettler Toledo, Germany), calibrated daily, was used.

For measurements of osmotic pressure a semi-micro osmometer (Knauer, Germany) was used. The semi-micro osmometer was calibrated before use over the range from 0 mOsmol/kg to 400 mOsmol/Kg with purified water and appropriate calibration solutions.

Methods

Measurements of the osmotic pressure (osmolality)

The measurement of the osmotic pressure of both dissolution media was performed according to the 4th European Pharmacopoeia (Ph. Eur. 4) monograph, which specifies the freezing point method. The osmotic pressure of each dissolution medium was determined in triplicate.

Buffer capacity

The buffer capacity was determined according to the titration method and reported as the volume of strong acid added to one liter buffer solution to lower its pH by one unit. This was determined by adding 1 N hydrochloric acid solution drop-wise from a 25 mL buret into the 100 mL sample buffer solution, with constant stirring, until the pH decreased by one pH unit. The investigations were performed in triplicate for each dissolution medium.

In vitro dissolution

In vitro dissolution tests were performed with the paddle apparatus in 500 mL dissolution media at 75 rpm. The dissolutions were run in accordance with the requirements of Ph. Eur. 4/ USP 26. Pre-measured volumes of dissolution media were degassed and then immediately transferred into the vessels. The tablets were placed in the vessels at the beginning of each test and the stopwatch was started simultaneously. Vessels were covered during the run with plastic covers in order to minimize evaporation. The temperature in the vessels was maintained at 37±0.5°C throughout each dissolution run. The revolution rate was checked during all tests. Samples were removed manually using 5 mL syringes fitted with stainless tubing to ensure reproducibility of the sampling location. Samples were filtered immediately after sampling through 45 µm PTFE filter sticks.

The parameters used for the *in vitro* dissolution testing are described in Table 3.

Analysis

The percent of drug dissolved at each sampling time point was determined for each drug product by UV analysis. The UV measurements were performed in all cases at the wavelength of the relative maximum absorption of the drug substance, which was determined prior to each analysis.

Calculation of % dissolved drug

Reference standard solutions were prepared by separately weighing portions of the chemical reference standard substance. Reference standard solutions were prepared and processed using the same method as for the preparation of the sample solutions. In order to avoid potential analytical interferences, standard solutions were prepared with aliquots of the same medium used for dissolution testing.

For each individual dissolution media UV-spectra between 200 nm and 400 nm of a typical blank solution, a typical sample solution, and the chemical reference standard solution were recorded. The relative maximum of the absorption spectrum was determined. By comparison of the recorded spectra, interference between the drug and the excipients dissolved in the medium was evaluated.

Concentrations of drug were calculated according to the external standard method with a linear calibration plot. The linearity of the calibration function was confirmed for each individual analytical sequence at three concentrations within the working range. The in-study validation for each sequence also included a check of the repeatability by six-fold measurement of reference standard solutions.

Results and Discussion

Osmotic pressure and buffer capacity

The results of the investigations on the osmotic pressure and buffer capacity of both dissolution media are given in Table 4.

The two media are practically identical with respect to pH, buffer capacity, ionic strength and osmotic pressure.

Table 4: Osmolarity, ionic strength and buffer capacity of the two buffers

Medium	Osmolarity [mOsmol/kg]	Ionic strength [mol/L]	Buffer capacity [mEq/L/pH unit]
Simulated Intestinal Fluid, pH 6.8 (SIFsp); USP 26	113	0.0720	18.4±0.2
Phosphate Standard Buffer pH 6.8 (IntPh 3)	115	0.0753	18.6±0.1

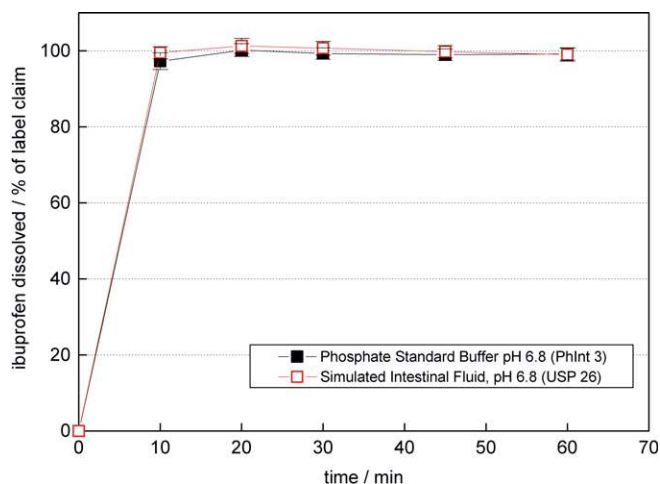


Figure 1.1. In vitro dissolution behavior of IBU-Product A film coated tablets. Arithmetic means (n=6) are shown, standard deviations were within the symbols.

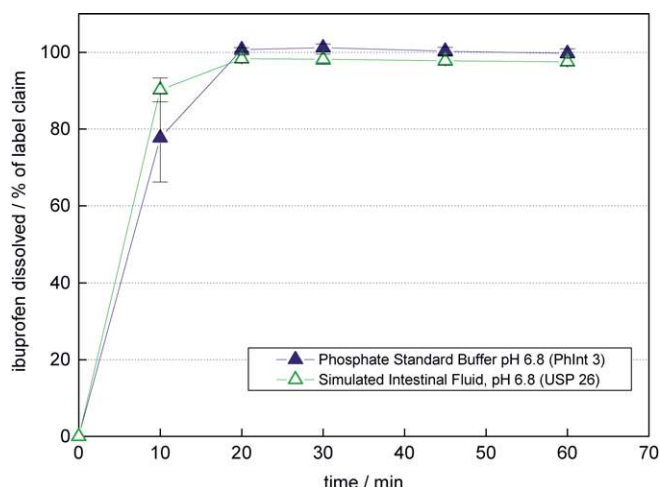


Figure 1.2. In vitro dissolution behavior of IBU-Product B sugar coated tablets. Arithmetic means (n=6) with standard deviations are shown.

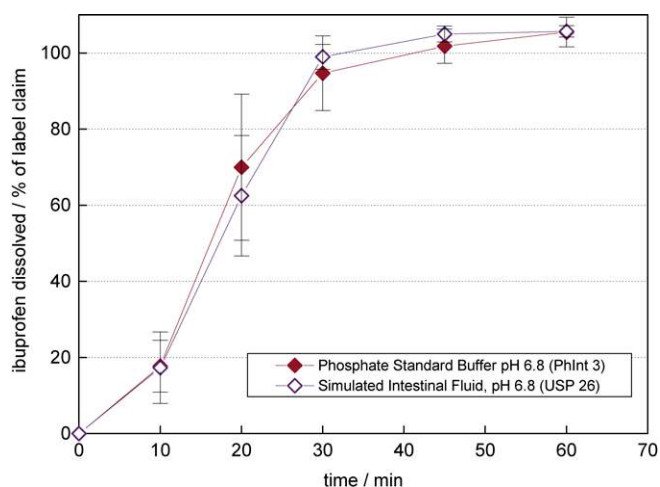


Figure 1.3. In vitro dissolution behavior of IBU-Product C capsules. Arithmetic means (n=6) with standard deviations are shown.

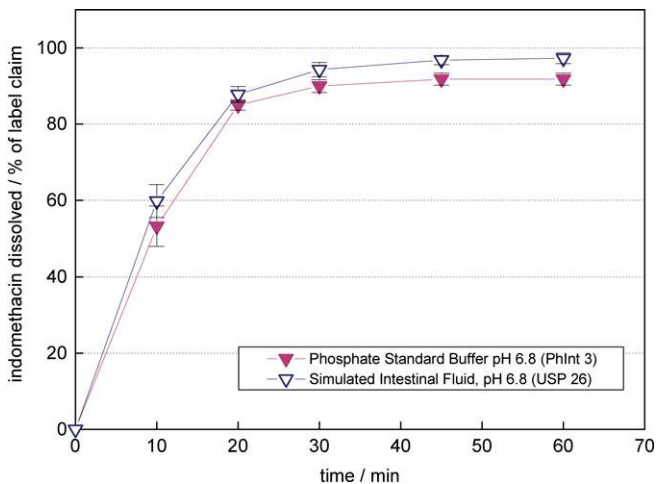


Figure 2.1. In vitro dissolution behavior of INDO-Product A uncoated tablets. Arithmetic means (n=6) with standard deviations are shown.

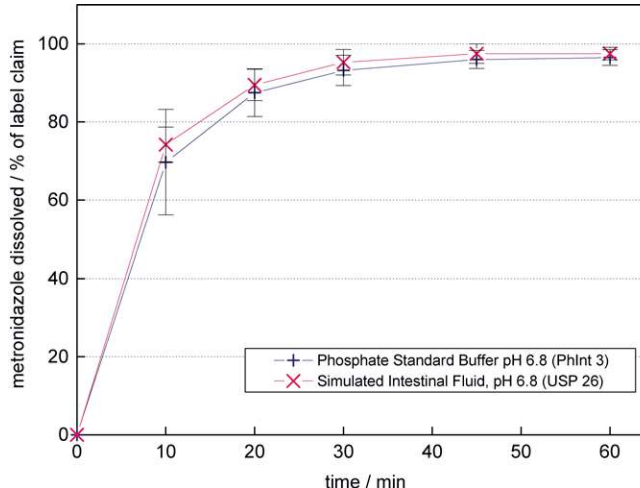


Figure 3.1. In vitro dissolution behavior of METRO-Product A film coated tablets. Arithmetic means (n=6) with standard deviations are shown.

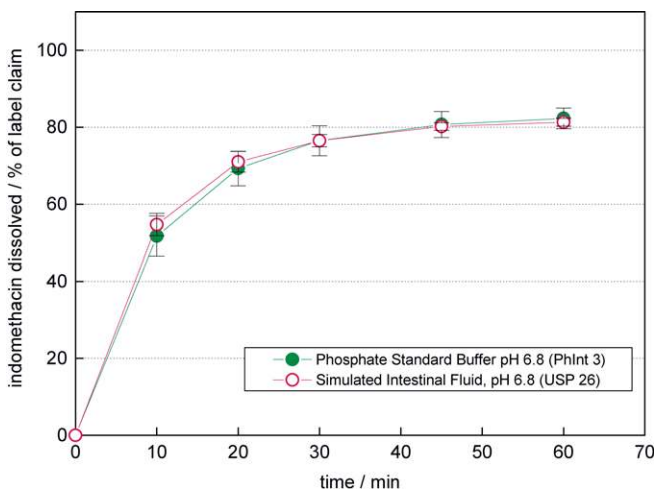


Figure 2.2. In vitro dissolution behavior of INDO-Product B capsules. Arithmetic means (n=6) with standard deviations are shown.

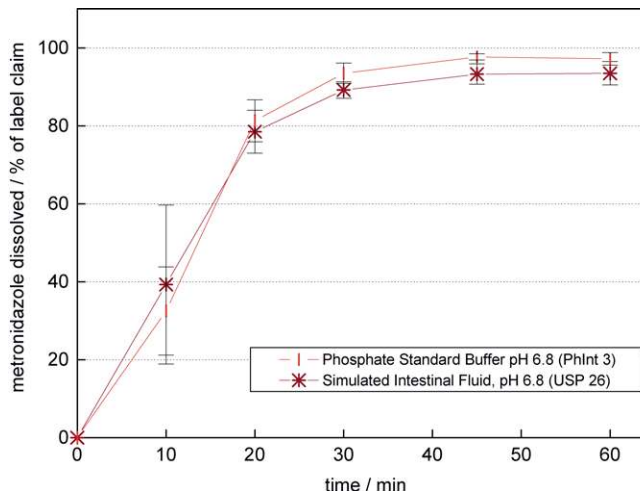


Figure 3.2. In vitro dissolution behavior of METRO-Product B uncoated tablets. Arithmetic means (n=6) with standard deviations are shown.

Therefore, dissolution behavior is expected to be similar in the two media for majority of products. Only in a specific case where the cation (Na^+ vs. K^+) may have an influence, would a difference in results be anticipated.

Dissolution Characterization

Excipient interference:

Interference of excipients did not exceed 5% for any of the products tested. Any interference was taken into consideration in calculating % dissolved, as described in the methods section.

Drug products containing ibuprofen

Mean *in vitro* dissolution behavior of various dosage forms containing ibuprofen under the two testing conditions are depicted in Figures 1.1 to 1.3.

No difference in the dissolution behavior of IBU-Product A film coated tablets was observed and dissolution was

complete in both media within 10 minutes. The dissolution profiles are superimposable and meet USP specifications.

A small difference in the initial dissolution behavior of IBU-Product B sugar coated tablets was observed. In both media the dissolution of IBU-Product B sugar coated tablets was complete within 20 minutes. The dissolution profiles are similar and both meet USP dissolution specifications.

Little or no difference in the dissolution behavior of IBU-Product C capsules was observed. Even though the dissolution profiles for this product are quite different from the other two products, they meet USP specifications. The difference is attributable to the dissolution time of the capsule shell.

Drug products containing indometacin

The *in vitro* dissolution behavior of two indometacin drug products (one uncoated tablet and one capsule

formulation) under the two testing conditions are depicted in Figures 2.1 and 2.2.

Little difference in the dissolution behavior of INDO-Product A uncoated tablets was observed. The amount of the drug dissolved in Phosphate Standard Buffer pH 6.8 was slightly lower in comparison to that in SIFsp, but the shape of the dissolution profiles are similar and both profiles meet USP specifications.

No difference in the dissolution behavior of INDO-Product B capsules was observed. The dissolution profiles are superimposable and meet USP specifications.

Drug products containing metronidazole

The *in vitro* dissolution behavior of two metronidazole drug products (one film coated and one uncoated tablet formulation) under the two testing conditions are depicted in Figures 3.1 and 3.2.

Little difference in the dissolution behavior of METRO-Product A film coated tablets in the two media was observed. The shapes of the dissolution profiles are similar and both profiles meet USP specifications.

Slight differences in the dissolution behavior of METRO-Product B uncoated tablets were observed. The amount of the drug dissolved at 10 minutes in Phosphate Standard Buffer pH 6.8 (IntPh 3) was slightly lower and after 20 minutes slightly higher in comparison to that in SIFsp. However, both profiles easily meet USP specifications.

Conclusions

Dissolution results were found to be similar in the three cases investigated. This is not surprising, since pH, buffer

capacity, ionic strength and osmolarity of the two buffers are almost identical. In fact, the distinguishing factor between the two media rests solely in the ratio of sodium to potassium. The results of these studies suggest that changing the proportions of these two cations in the medium makes little or no difference to the dissolution results. For the purposes of dissolution testing of conventionally formulated immediate release solid oral dosage forms, the Phosphate Standard Buffer pH 6.8 of the IntPh 3 can be generally considered interchangeable with SIFsp of the USP 26. Only in specific cases where solubility is known to be affected by the cation does one need to be wary about substitution of the two cations for each other. In an upcoming report, we will address this issue using the dissolution of phenytoin sodium capsules in sodium lauryl sulfate containing media as a case example.

References

1. United States Pharmacopeia and National Formulary, United States Pharmacopeial Convention Inc., Rockville, MD, USA, 25th Edition (2002) and 26th Edition (2003)
2. The International Pharmacopoeia, WHO, Geneva, Switzerland, 3rd Edition, Volume 1 (1979) and Volume 5 (2003)
3. European Pharmacopoeia, European Directorate for the Quality of Medicines of the Council of Europe (EDQM), Strassbourg, France, 4th Edition (2001)