

*Conceptual approach to data integration and modelling, implemented using the megNet® software. Both statistical and semantic models are utilized to enable systemic integration of data across multiple levels. The platform also enables integration of models and knowledge across multiple species.*

tion AND females", from which a network of relations in clinical and biological databases is built and visualized.

#### Perspectives

Our conceptual space strategy, implemented using the megNet® environment, has already demonstrated its potential in clinical applications. We believe that our approach will be very useful in building complex in silico models at the level of human physiology, making mappings across multiple levels of biological organization and

across multiple knowledge domains a feasible task.

One of the emergent challenges in life science and medical knowledge management is how to deal with the dynamics in biological systems, that is, how to encode the inherent dynamic properties of biological systems for the purpose of data mining and modelling. It is obvious that modelling at all levels, from quantum processes to physiology and environment, is computationally unfeasible. We believe that the conceptual spaces

approach could help in establishing the relevant components of the system to be included in the models.

Since conceptual spaces are a powerful approach to build metaphors across different knowledge domains, one could also envision the applications of the approach outside the life science domain. We have recently initiated one such project, aiming to use agent-based modelling of biological cells in order to develop more flexible computing tools.

#### Links:

Quantitative Biology and Bioinformatics group at VTT: <http://sysbio.vtt.fi/>  
 VISUBIOMED project: <http://sysbio.vtt.fi/visubiomed/>  
 TRANSCENDO project: <http://sysbio.vtt.fi/transcendo/>  
 SYSDIPP project: <http://sysbio.vtt.fi/sysdipp/>  
 Helsinki University of Technology, Laboratory of Computer and Information Science: <http://www.cis.hut.fi/>

#### Please contact:

Matej Orešič  
 VTT Technical Research Centre of Finland  
 Tel: +358 20 722 4491  
 E-mail: [matej.oresic@vtt.fi](mailto:matej.oresic@vtt.fi)

## In Vivo Microscopy for Real-Time Structural and Functional Cellular Imaging

by Tom Vercauteren, Aymeric Perchant and Nicholas Ayache

*Fibered confocal microscopy allows the acquisition of in vivo and in situ images at the cellular level, in combination with standard endoscopic procedures or needle biopsies for solid organs. This makes it a promising tool for clinical molecular imaging, an activity aiming at in vivo characterization and measurement of biological processes. Confocal microscopy images represent a new source of information for developing patient-specific digital models that integrate knowledge of cellular dynamics. This is also an unrivaled technique for refining digital patient models down to the microscopic level.*

Current medical practice is moving towards more quantitative and personalized methods of diagnosis and decision-making. This has led to a need for more complex and detailed patient-specific digital models. The resolution at which the model is described is important and may vary from macroscopic to microscopic and even molecular scales, ideally through multiscale descriptions. As these individualized models start to inte-


grate knowledge of cellular dynamics, it becomes crucial to acquire patient-specific information at the cellular level.

The goal of the cooperation between the Asclepios research group at INRIA Sophia Antipolis and the company Mauna Kea Technologies (MKT), Paris, is to develop advanced image analysis tools capable of extracting the

pertinent information from in vivo confocal microscopy images. This is a first step towards the integration of individualized cellular dynamics information in the physiological models we develop at Asclepios.

#### Fibered Confocal Microscopy

Mauna Kea Technologies' Cellvizio® family of endo-microscopes rely on a unique approach for the in vivo and in

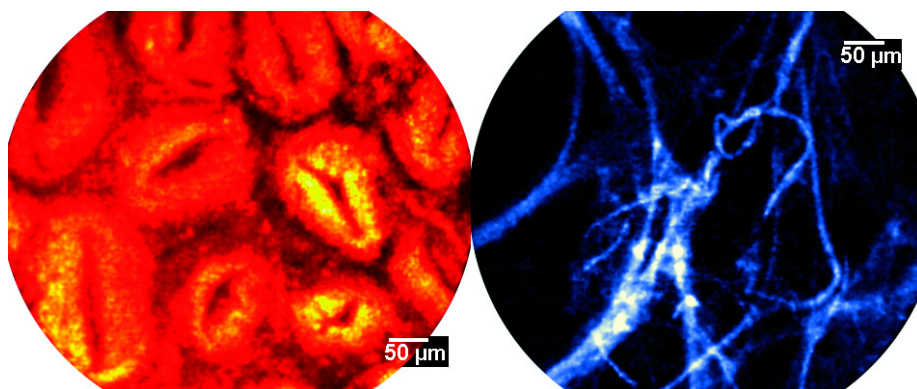

*Figure 1: Fibered confocal microscopy allows the clinician to perform micron-scale resolution optical biopsies and to diagnose the state of a tissue in vivo and in real time. Left: Normal gastric pits. Courtesy: PD. Dr. A. Meining, Klinikum rechts der Isar, Munich. Right: Alveolar Network. Courtesy: Pr. L. Thiberville, CHU Charles Nicolle, Rouen.*


situ exploration of living organisms. Simply stated, this involves putting a microscope objective at the end of an ultra-thin, three meter long optic fiber.

The Cellvizio® acquires image sequences of microscopic resolution, displays them in real time and enables live measurements. The confocal miniprobe, one of the Cellvizio®'s three components, is a custom-built highly advanced optical imaging probe incorporating a proprietary fiber bundle and objective lens technology. Coupled to the laser scanning unit, the sophisticated image processing software renders real-time dynamic image sequences with a lateral resolution as fine as 2.5µm and a field of view of up to 600µm at 12 frames per second.

This specific imaging modality raises certain image processing problems. For example the raw data generated by the instrument shows geometric distortions and a non-uniform honeycomb pattern due to the modulation of the fiber optic bundle. This makes the data impracticable for automated analysis if left untreated. Algorithms that take on the image reconstruction task in real time have thus been developed in order to provide high-quality, smooth-motion video sequences.

### Bridging the Gap between Micro- and Macroscopic Scales

Fibered confocal microscopy can unveil in real time the cellular structure of the observed tissue. However, as interesting as dynamic sequences may be during the time of the medical procedure or biological experiment, it is necessary to have an efficient representation of the entire imaged region. Mosaicing techniques can be used to provide this representation. However classical algo-

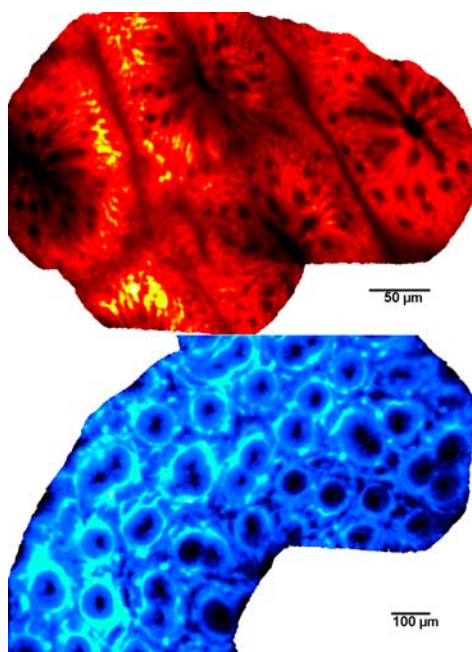

*Figure 2: Mosaic reconstructions can increase the field of view while enhancing the image definition. Top: Healthy colonic mucosa (26 images). Bottom: Microscopic colitis (98 images). Courtesy: PD. Dr. A. Meining, Klinikum rechts der Isar, Munich.*

rithms do not take into account the characteristics of fibered confocal microscopy, namely motion distortions, irregularly sampled frames and non-rigid deformations of the imaged tissue. This is why specific novel algorithms were developed.

With these mosaics, quantitative and statistical analysis becomes possible over a wide field of view. The quality of the data that can be used to build cellular-level patient-specific models is thus improved. Moreover, mosaicing for microscopic images is a means of bridging the gap between microscopic and macroscopic scales. The combination of Cellvizio® and image processing therefore provides a unique way to build patient-specific digital models at the microscopic scale.

### Links:

<http://www.maunakeatech.com>  
<http://www-sop.inria.fr/asclepios>

### Please contact:

Tom Vercauteren, Aymeric Perchant  
 Mauna Kea Technologies, Paris, France  
 Tel: +33 1 48 24 11 43  
 E-mail: [tom@maunakeatech.com](mailto:tom@maunakeatech.com),  
[aymeric@maunakeatech.com](mailto:aymeric@maunakeatech.com)

Nicholas Ayache  
 Asclepios project team, INRIA, Sophia Antipolis, France  
 Tel: +33 4 92 38 76 60  
 E-mail: [nicholas.ayache@sophia.inria.fr](mailto:nicholas.ayache@sophia.inria.fr)